

PCT


WORLD INTELLECTUAL PROPERTY ORGANIZATION  
International Bureau



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification 6 : <b>G06F 15/00, A61B 5/00, G01N 33/48, 31/00, 33/00, G06K 9/00</b>		<b>A1</b>	(11) International Publication Number: <b>WO 95/06909</b>
(21) International Application Number: <b>PCT/US94/10001</b>		(43) International Publication Date: <b>9 March 1995 (09.03.95)</b>	
(22) International Filing Date: <b>31 August 1994 (31.08.94)</b>		<p>(81) Designated States: AU, CA, European patent (AT, BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE).</p> <p><b>Published</b>  <i>With international search report.          Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.</i></p>	
(30) Priority Data: <b>08/116,131      2 September 1993 (02.09.93)      US</b>			
(71) Applicant: COULTER CORPORATION [US/US]; 11800 S.W. 147 Avenue, Miami, FL 33196-2500 (US).			
(72) Inventors: DIAMOND, Lawrence, W.; Uni-Center, Apartment #619, Luxemburger Strasse 124, D-50939 Cologne (DE). NGUYEN, Doyen, T.; Uni-Center, Apartment #619, Luxemburger Strasse 124, D-50939 Cologne (DE).			
(74) Agent: WINBURN, John, T.; Coulter Corporation, Patent Dept., M/C 32-A02, P.O. Box 169015, Miami, FL 33116-9015 (US).			

(54) Title: DECISION SUPPORT SYSTEM AND METHOD FOR DIAGNOSIS CONSULTATION IN LABORATORY HEMATOPATHOLOGY

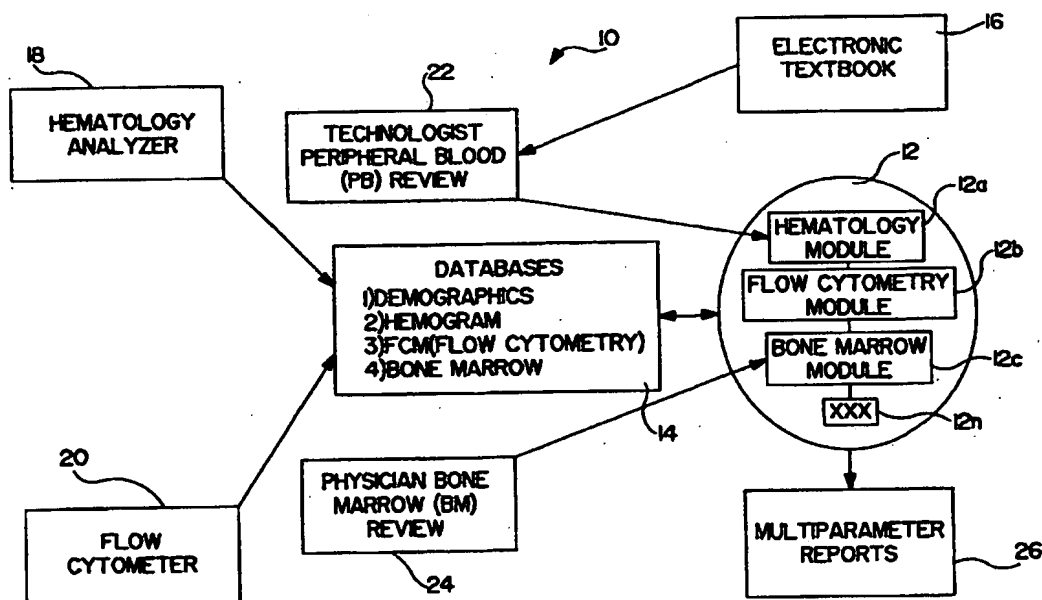

(57) Abstract

A system (10) and method of diagnosis provides that data from a hematology analyzer, e.g., hematology analyzer (18) and/or other blood testing instrument or test, is input and compared with patterns corresponding to specific patient conditions. The matched patterns are arranged in a hierarchy in accordance with predetermined rules. Initially the patterns matched are displayed, optionally with recommendations for additional testing. If additional testing is done, a comparison of the original test data and new data is done, the patterns matched and at least one highest ranked pattern is identified, optionally with a differential diagnosis.

**FOR THE PURPOSES OF INFORMATION ONLY**

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	GB	United Kingdom	MR	Mauritania
AU	Australia	GE	Georgia	MW	Malawi
BB	Barbados	GN	Guinea	NE	Niger
BE	Belgium	GR	Greece	NL	Netherlands
BF	Burkina Faso	HU	Hungary	NO	Norway
BG	Bulgaria	IE	Ireland	NZ	New Zealand
BJ	Benin	IT	Italy	PL	Poland
BR	Brazil	JP	Japan	PT	Portugal
BY	Belarus	KE	Kenya	RO	Romania
CA	Canada	KG	Kyrgyzstan	RU	Russian Federation
CF	Central African Republic	KP	Democratic People's Republic of Korea	SD	Sudan
CG	Congo	KR	Republic of Korea	SE	Sweden
CH	Switzerland	KZ	Kazakhstan	SI	Slovenia
CI	Côte d'Ivoire	LI	Liechtenstein	SK	Slovakia
CM	Cameroon	LU	Luxembourg	SN	Senegal
CN	China	LV	Latvia	TD	Chad
CS	Czechoslovakia	MC	Monaco	TG	Togo
CZ	Czech Republic	MD	Republic of Moldova	TJ	Tajikistan
DE	Germany	MG	Madagascar	TT	Trinidad and Tobago
DK	Denmark	ML	Mali	UA	Ukraine
ES	Spain	MN	Mongolia	US	United States of America
FI	Finland			UZ	Uzbekistan
FR	France			VN	Viet Nam
GA	Gabon				

**DECISION SUPPORT SYSTEM AND METHOD FOR DIAGNOSIS  
CONSULTATION IN LABORATORY HEMATOPATHOLOGY**

**Technical Field**

This invention relates to a system and method for conducting laboratory hematopathology diagnosis. By "hematopathology" is meant the medical science  
5 concerned with diseases of the blood and related tissues, including but not limited to spinal fluid, lymph nodes and bone marrow. Hematopathologists are generally considered to include such persons as hematologists, pathologists, immunologists and  
10 related scientists.

More specifically, the invention relates to a system and method for conducting hematopathology diagnosis based on results output from a hematology analyzer, flow cytometry (FCM) analyzer, and review  
15 of microscopic slides of bone marrow preparations. The invention provides a patterned approach to peripheral blood, FCM, and bone marrow analysis and serves to guide technologists and technicians through the complexities of hemogram interpretation,  
20 morphologic observations and complementary laboratory studies. Furthermore, the invention guides physicians in the interpretation of FCM studies and provides a means of generating multiparameter interpretive reports of hematology  
25 specimens. A logical mechanism of diagnostic reasoning is implemented through the use of a computer.

-2-

Background Art

In medicine, diagnostic reasoning is the process of converting data obtained from examining a patient into names of diseases, which are of themselves conceptual medical entities that identify or explain the abnormalities in the observed evidence. In practice in a clinical setting, diagnostic reasoning often does not proceed in a direct continuous manner towards a precise disease identification, but is often interrupted by the performance of laboratory tests. Diagnostic consultation thus requires a strategy for evaluating evidence, recognizing patterns, ordering adjunctive tests and ultimately choosing the names of diagnostic entities.

A common screening test in clinical medicine is the hemogram or "complete blood count" (CBC). Analysis of hemogram data generated by sophisticated, modern hematology analyzers capable of automated "white blood cell" (WBC) differentials is routinely performed for most hematologic as well as non-hematologic disorders. Correct interpretation of the data will then indicate whether or not examination of the peripheral blood smear is necessary, and the hemogram findings together with the impression formed from scanning the blood smear, indicate when a manual WBC

-3-

differential should be performed. Interpretation of the hemogram and peripheral blood smear findings leads to what is referred to conventionally as a diagnostic "station," an intermediate hypothesis during the diagnostic reasoning process. This "station" or "pattern" as referred to hereinafter, results in a narrow group of differential diagnoses and should provide an indication of which additional tests will confirm the diagnosis with minimum costs, delay and discomfort to the patient. Presently, conventional hematology analyzers, such as the commercially available Coulter® STKS hematology analyzer sold by Coulter Corporation of Miami, Florida, provide cell size distribution curves for red blood cells, WBC and platelets (PLT). The Coulter® STKS device analyzes the WBC populations by a combination of volume, conductivity and light scatter (VCS) measurements. An automated five-part WBC differential is then performed and suspect flags are generated for potential morphologic abnormalities by means of cluster analysis.

Typically, the correct interpretation of hematologic patterns, as well as the morphologic findings and the results of such tests as immunologic analysis of leukocytes by flow cytometry, is often difficult for laboratory personnel and physicians lacking specific training

-4-

and expertise in the fields of hematology and hematopathology.

To facilitate such diagnosis, a number of computerized decision support programs have been developed to aid in the interpretation of hemogram findings. Of particular interest is that such programs have been developed with reference to the diagnosis of anemia. Often the approach employed by these systems include Bayesian probability as disclosed in "Decision-making systems (DMS) applied to hematology. diagnosis of 180 cases of anemia secondary to a variety of hematologic disorders" by Sultan C., Imbert M., Priolet G.; Hematol Pathol 1988; 2:221-8; statistical pattern recognition as disclosed in "The use of an expert system in the clinical laboratory as an aid in the diagnosis of anemia" by Blomberg D.J., Ladley J.L., Fattu J.M., Patrick E.A.; AM J Clin Pathol 1987; 87:608-13; multi-variate statistical analysis as disclosed in "Evaluation of BCDE, a microcomputer program to analyze automated blood counts and differentials" by Bates J.E., Bessman J.D.; AM J Clin Pathol 1987; 88:314-23, and rule-based production systems as disclosed in "An expert system for the interpretation of full blood counts and blood smears in a hematology laboratory" by Tolmie C.J., du

-5-

Plessis J.P., Badenhorst P.N.; Artif Intel Med 1991;  
3:271-85.

When such systems have been evaluated with a limited number of anemia cases, performance was  
5 judged satisfactory in about 84%-89% of the cases. While encouraging, the success numbers in terms of arriving at a correct diagnosis falls short of what is desired in an actual diagnostic environment.

Thus, in accordance with the invention, the  
10 problems of the prior art are avoided by employing a system and method which does not rely upon probability/statistics methodology to arrive at a differential diagnosis or interpretation of test results. Instead, a hierarchical classification of  
15 blood, immunophenotype and bone marrow patterns is established and a heuristic algorithm is applied to eliminate patterns based on data obtained from hematopathology analysis. For purposes of this disclosure, by "differential diagnosis" is meant a  
20 list of possible diagnoses or diagnostic possibilities. Further, by "interpretation of test results" is meant a diagnostic impression based on available data.

#### Disclosure of the Invention

25 A method of interpreting hematopathology test results to arrive at a diagnosis of a patient

-6-

characterized by retrieving data from at least one  
hematopathology test for a patient and storing said  
hematopathology test data in a database; comparing  
said hematopathology test data with at least one set  
5 of a predetermined number of hematopathology  
patterns, each of said hematopathology patterns  
corresponding to different patient conditions and  
said patterns arranged in a predetermined hierarchy,  
to determine which patterns match said  
10 hematopathology test data in said database; and  
displaying the patterns matched, with  
recommendations corresponding to the matched  
patterns for one of no additional testing or one of  
a plurality of identified additional testing, the  
15 results of which are likely to lead to a diagnosis  
of the patient's condition.

The method of interpreting hematopathology  
test results to arrive at a diagnosis of a patient;  
characterized by retrieving data from at least one  
20 hematopathology test for a patient and storing said  
hematopathology test data in a database; comparing  
said hematopathology test data with at least one set  
of a predetermined number of hematopathology  
patterns, each of said hematopathology patterns  
25 corresponding to different patient conditions and  
said patterns arranged in a predetermined hierarchy  
to determine which patterns match said



-7-

hematopathology test data in said database; and identifying at least one highest ranked pattern in said hierarchy to arrive at a differential diagnosis of the patient's condition.

5           A system for interpreting hematopathology test results to arrive at a diagnosis of a patient, characterized by database storage means for storing a database of test results from hematopathology tests for specific patients; knowledge  
10 base/inference engine means for storing at least one set of a predetermined number of patterns corresponding to specific patient conditions, including a normal condition, with said patterns arranged in a predetermined hierarchy, and for  
15 storing instructions for comparing said patterns to said database of test results for determining which patterns are matched; comparing means for comparing data from said tests for a specific patient with said patterns to determine which patterns are  
20 matched and which pattern is highest in the hierarchy; input means for inputting said data from test results into said database; and output means for outputting the results of said comparison by said comparing means, said results being comprised  
25 of at least one of the highest ranked pattern, a differential diagnosis and recommendations for additional testing.

-8-

In one aspect, the invention is directed to a method of interpreting hematopathology test results to arrive at a diagnosis of a patient. As noted previously, hematopathology tests can include  
5 hematology, lymph node, bone marrow, flow cytometry, spinal fluid and other types of testing generally done to diagnose diseases of the blood and related tissues.

In accordance with its broadest aspect, the  
10 method includes retrieving data from at least one hematopathology test for a patient and storing the data in a database. The hematopathology test data is then compared to at least one set of a predetermined number of hematopathology patterns.  
15 Each of the hematopathology patterns corresponds to different patient conditions and they are arranged in a predetermined hierarchy. A comparison is done to determine which patterns match the hematopathology test data in the database. The  
20 matched patterns are thereafter displayed with recommendations corresponding to the matched patterns including either no additional testing (for example, if the patient condition is indicated to be normal) or a recommendation for one of a plurality  
25 of additional testing, the results of which are likely to lead to a diagnosis of the patient's condition.

-9-

As further refined, the method includes conducting any additional testing recommended, retrieving the data therefrom and storing it in the database with the hematopathology test data. A  
5 second comparison is then done with the original as well as the new data with the patterns, to determine which patterns match the test data in the database. At least one highest ranked pattern is identified to arrive at a differential diagnosis of the patient's  
10 condition.

As may be appreciated, in the ordinary course of diagnosis, the first hematopathology test conducted is a hematology analysis which may include a visual examination of a blood smear. This  
15 comparison may yield fairly conclusive results, may also lead to a suggestion of additional testing for a more definitive differential diagnosis or interpretation of the patient's test results. Thus, the additional testing recommended from the  
20 hematology test may include, but is not limited to flow cytometry or bone marrow testing. The flow cytometry or bone marrow data would then also be entered into the database and a follow-on comparison made of all the data to patterns, in the case of  
25 flow cytometry to a set of flow cytometry patterns, and in the case of bone marrow testing, to a set of bone marrow patterns, which would include as part of

-10-

their identification any previous differential diagnosis resulting from the hematology testing.

Although hematology testing is indicated as typically being done first, it will be appreciated  
5 that either type of testing can be done as a stand-alone procedure, with the flow cytometry and bone marrow testing being more specific in results than the hematology testing. Although from the flow  
10 cytometry or bone marrow test results it is possible to arrive at a fairly specific diagnosis, it will also be generally recommended that a hematology test be done before reaching a conclusive diagnosis. In the case of flow cytometry it will be appreciated the multiple patterns each of equal and highest rank  
15 may be identified. Thus, this is why, for example, a hematology test may be indicated.

In another aspect, the invention is also directed to a system for interpreting  
hematopathology test results to arrive at a  
20 diagnosis of a patient. The system includes database storage for storing a database of test results from hematopathology test for specific patients. A knowledge base stores at least one set of a predetermined number of patterns corresponding  
25 to specific patient conditions, including at least a normal condition, and with the patterns arranged in a predetermined hierarchy. The knowledge base

-11-

and an inference engine contain instructions for comparing the patterns to the database of test results for determining which patterns are matched. In the system of the invention, the knowledge base and inference engine are linked in the programming as implemented. A comparator serves to conduct the comparison from the tests for a specific patient with the patterns to determine which patterns are matched and which pattern is highest in the hierarchy. An input device serves to input data from the test results into the database and an output device serves to output the results of the comparison.

In accordance with a more specific aspect of the invention, there is provided a method of interpreting hematology results to arrive at a differential diagnosis for a patient with recommendations for additional diagnostic tests in the case of peripheral blood analysis. The method includes the steps of establishing a knowledge base of a predetermined number of exhaustive but not mutually exclusive hematology patterns. By the term "exhaustive number of patterns" is meant all of the most commonly occurring patterns, including a normal pattern, and one pattern which corresponds to abnormal conditions but not specific to any condition. The nonspecific pattern accounts for no

-12-

more than about 5% of abnormal conditions showing what are typically nonsignificant abnormalities for diagnostic purposes. Further, by the term "knowledge base" is meant a set of defined patterns and programmed instructions for matching the patterns to input data for constructing a report from the matching.

One of the patterns in the knowledge base corresponds to a normal patient condition and each of the remaining patterns except one is representative of a corresponding differential diagnosis of specific conditions. The last predefined pattern in the peripheral blood module is designated "non-specific" and, as previously noted, includes cases which do not meet the defined criteria for any of the other aforementioned patterns. Data is retrieved from at least a first hematology test for a specific patient. The data retrieved is then compared with the predetermined number of hematology patterns in the knowledge base to determine which hematology patterns match the data.

An input device is provided for inputting data from hematology tests for a specific patient. A comparing device serves to compare the data input by the input device with the predetermined number of hematology patterns in the database. An arranging

-13-

or sorting device serves to arrange the hematology patterns which match the input data in a predetermined hierarchy. Since the patterns are not mutually exclusive, a predetermined hierarchy of the patterns is established.

Diagnostic reasoning proceeds according to the pattern matched which is highest in the hierarchy. A display displays the pattern highest in the hierarchy for indicating a differential diagnosis of the patient's condition. Computer program modules are implemented in the form of what are conventionally known as "if...then" statements to generate recommendations for additional tests which may help to establish the final diagnosis of the patient's condition. These recommendations are displayed along with the peripheral blood pattern and the differential diagnosis.

In another aspect, the invention relates to a system for interpreting FCM results for generating an interpretation with a differential diagnosis in the case of FCM analysis. The system optionally includes database storage for both the results of the peripheral blood findings and the peripheral blood interpretation described above. The system also includes a knowledge base of a predetermined number of exhaustive but not mutually exclusive FCM

-14-

patterns. Each pattern except one corresponds to a specific cell type (lineage). One pattern, designated "non-specific" corresponds to conditions which do not meet the definitions of the specific patterns because one or more FCM markers was not included in the panel of markers studied in a given patient specimen.

An input device is provided for inputting data from FCM tests for a specific patient. A comparing device serves to compare the data input by the input device with the predetermined number of hematology patterns in the knowledge base. An arranging and sorting device serves to arrange the FCM patterns which match the input data. Database programmed functions are provided to look up any corresponding peripheral blood information stored in the database on the patient. Computer program modules are implemented in the form of "if...then" statements to use the input data and the data retrieved from the database in the interpretation of the findings. A display displays the FCM pattern(s) matched for indicating a differential diagnosis of the patient's condition and the remaining parts of the interpretation as generated by the computer program modules.

In another aspect, the invention relates to a system for aiding in the interpretation of bone



-15-

marrow microscopic slides and generating a multiparameter report of the findings. As described above, the system optionally includes database storage for the peripheral blood findings and interpretation. The system also optionally includes database storage for the FCM findings and interpretation as well as the bone marrow findings and interpretation. In addition, the system also includes a knowledge base with a predetermined number of exhaustive but not mutually exclusive bone marrow patterns. Each pattern except one is representative of a corresponding differential diagnosis of specific conditions. One pattern, designated "non-specific" corresponds to conditions which do not meet the definition of the specific patterns.

An input device is provided for inputting data from observation of bone marrow microscopic slides by a physician (and a trained technologist, where appropriate). A comparing device serves to compare the data input by the input means with the predetermined number of bone marrow patterns in the knowledge base. An arranging and sorting device serves to arrange the bone marrow patterns which match the input data.

Since the patterns are not mutually exclusive, a predetermined hierarchy of the patterns

-16-

is established. Diagnostic reasoning proceeds according to the pattern matched which is highest in the hierarchy. Database programmed functions are provided to look up any corresponding peripheral blood and FCM information stored on the patient. Computer program modules are implemented in the form of "if...then" statements to use the input data and the information retrieved from the database in the interpretation of the bone marrow findings. A display displays the appropriate peripheral blood information retrieved from the database, the highest bone marrow pattern matched in the hierarchy, the bone marrow morphologic findings input by the user, the FCM interpretation retrieved from the database (if FCM was performed on the a peripheral blood or bone marrow specimen from this patient) and the final diagnosis determined by the computer program modules.

In a more specific aspect, the system is implemented on a computer including a keyboard and a "direct access storage device" (DASD) connected thereto, such as a hard disk drive. In addition, the input devices can include an electrical interface directly connected to the computer and adapted for direct connection to a hematology analyzer, and optionally to a flow cytometry analyzer, for having the results of blood cell

-17-

analysis, and flow cytometry analysis, conducted thereby input directly into the computer. The comparing devices are implemented in conventional "comparator/sorter" hardware components with  
5 associated instructions in such a computer.

#### Brief Description of Drawings

The preferred embodiments of this invention will now be described by way of example, with reference to the drawings accompanying this  
10 specification in which:

FIG. 1 is a block diagram illustrating the general system in accordance with the invention, including the interaction with a technician and/or a physician operating the system for conducting  
15 diagnosis;

FIG. 2 is a block diagram of the typical architecture of a decision support system in accordance with the invention, including a knowledge base, inference engine, and a user interface;

20 FIG. 3 is a general block diagram illustrating a typical computer system on which the invention can be implemented;

FIG. 4 is a block diagram of a typical computer of the type shown in FIG. 3, illustrating  
25 the various components thereof;

-18-

FIGS. 5, 6 and 7 are reproductions of actual screens from the hematology module to illustrate the user interface and the ability of the technologist to overrule the system recommendations based on  
5 observation of the peripheral blood smear at scanning magnification.

FIGS. 8, 9, 10, 11 and 12 are reproductions of actual screens from the hematology module to illustrate the operation of this module in a patient  
10 with acute leukemia.

FIGS. 13, 14 and 15 are screen displays illustrating the operation of the flow cytometry, i.e., FCM module of the invention.

FIGS. 16, 17, 18 and 19 are screen displays  
15 illustrating the operation of the bone marrow module of the invention in the case of a patient with acute myeloid leukemia.

FIGS. 20 and 21 illustrate actual screen displays from the electronic textbook associated  
20 with the hematology module, also described as "Professor Petrushka."

FIG. 22 is a block diagram illustrating the steps taken in conducting the method in accordance with the invention.

-19-

Modes for Carrying out the Invention

Initially, it is noted that throughout this specification, the terms "hematology" and "peripheral blood" findings are used interchangeably and mean the same thing. Further, the "hematology module" is interchangeably referred to as the "peripheral blood module," or as shown in the screen displays in the figures, as "Professor Petrushka" or "Petrushka." The flow cytometry or FCM module is interchangeably referred to as "Professor Fidelio" or "Fidelio," and the bone marrow module is interchangeably referred to as "Belmonte" or "Professor Belmonte."

The decision support system in accordance with the invention is generally shown in FIG. 1 and designated by the numeral 10. Included in the system 10 are a number of modules 12 made up of programmed instructions, as will become more clearly evident from the following discussion, which operate within a computer to process data to arrive at a differential diagnosis based on the results of hematopathology analysis.

In a preferred aspect, the invention as presently contemplated includes three modules, a

-20-

hematology module 12a, a flow cytometry module 12b and a bone marrow module 12c for processing the results of different tests on a patient. The modules 12a, 12b and 12c interact with databases 14  
5 corresponding to each module and each module contains a knowledge base 13, as shown in FIG. 2, which illustrates the typical architecture of the decision support system of the invention, which includes a number of patterns each corresponding to  
10 a list of differential diagnoses.

In its least complex implementation, the system is implemented only with the hematology module 12a, the appropriate database 14 and the corresponding knowledge base 13, inference engine  
15 15, and a user interface 17, as shown in FIG. 2, which for interpretation of the results of hematology analyzer 18 as described hereafter.

For hematology analysis the patterns shown in the following Table I are stored in the knowledge  
20 base 13:

**TABLE I**  
**PERIPHERAL BLOOD PATTERNS**

1. Normal	10. Monocytosis
2. Reactive Neutrophilia	11. Bi-Pancytopenia
3. Leukoerythroblastic	12. Leukopenia
4. Abnormal Mononuclear Cell	13. Hypochromic and/or Microcytic
5. Erythrocytosis	14. Macrocytic
6. Thrombocytosis	15. Normochromic, Normocytic Anemia
7. Absolute Lymphocytosis	16. Anemia with Anisocytosis
8. Thrombocytopenia	17. Leukemoid
9. Eosinophilia	18. Non-specific

As shown in FIG. 3, a computer 30 of the system 10 is typically connected through a hardware interface 36 to a hematology analyzer 18. As also shown in FIG. 3, it is possible to implement the invention with more than a hematology analyzer 18 in a complete hematopathology diagnostic system. For example, it is also possible to connect a flow cytometry analyzer 20 through interface 36a to

-22-

computer 30. In this embodiment the flow cytometry module (FCM) will contain the corresponding knowledge base 13, inference engine 15, and user interface 17 as shown in FIG. 2, and the database 14 of FIG. 1 will contain a table of data appropriate for FCM data as shown in Appendix 3, and the knowledge base 13 contains flow cytometry patterns as shown in Table II below. Similarly, the system can also include a bone marrow module 12c, with the bone marrow module 12c patterns employed in a similar manner as discussed hereafter and shown in Table III following Table II below.

TABLE II

## FLOW CYTOMETRY IMMUNOPHENOTYPING PATTERNS

- |    |                    |                    |
|----|--------------------|--------------------|
| 15 | 1. Stem Cell       | 6. Immature T-cell |
|    | 2. Myeloid         | 7. Mature B-cell   |
|    | 3. Erythroid       | 8. Mature T-cell   |
|    | 4. Megakaryocytic  | 9. Plasma cell     |
|    | 5. Immature B-cell | 10. Non-specific   |

TABLE III

## BONE MARROW PATTERNS

- |    |  |                                  |
|----|--|----------------------------------|
| 20 | 1. Normocellular                                       | 7. Mononuclear cell infiltration |
|    | 2. Hypocellular  | 8. Foreign cell infiltration     |
|    | 3. Erythroid predominance/hyperplasia                  | 9. Granulomatous                 |
| 25 | 4. Hypercellular, myeloid                              | 10. Bi-hyperplasia               |
|    | 5. Megakaryocytic hyperplasia predominance/hyperplasia | 11. Non-specific                 |
|    | 6. Lymphocytosis/plasmacytosis                         |                                  |



-23-

Typically, the system 10, as shown in FIG. 3, will consist of a computer 30 connected to an input device such as a keyboard 34, and an output device such as a display 32. The computer 30 can also be  
5 directly connected through the electrical interface 36 to a hematology analyzer 18 through which the results of a test conducted therein can be directly input into the computer. Such a hematology analyzer 18 can be, for example, a Coulter® STKS hematology  
10 analyzer, such as sold by Coulter Corporation of Miami, Florida, which is presently commercially available and well known to those of ordinary skill in the art.

The specifics of the computer 30 on which the  
15 system 10 is implemented and which is conventional, are more clearly shown in FIG. 4 which shows the computer 30 including a conventionally arranged central processing unit 40 including "read only memory" (ROM) 44 which includes, for example, the  
20 necessary basic input/output system instructions such as employed in an IBM® personal computer under the name BIOS. "random access memory" (RAM) 42 serves to store data and instruction from a keyboard

-24-

or other input device, as well as from a "direct access storage device" (DASD) 46, such as a hard disk drive, or a floppy disk drive 48.

The input from the analyzer 18 will include  
5 WBC, RBC and platelet (PLT) population flags and suspect flags for morphologic abnormalities which should be looked for on the peripheral blood smear.

One suspect flag is for blasts which are known to those of ordinary skill in the art as being  
10 the most immature white blood cells. Blasts are found circulating in the blood in large numbers only in serious malignant conditions such as acute leukemia and related disorders. Therefore, a suspect blast flag means that the corresponding  
15 blood smear should be carefully reviewed. Other suspect flags are for: 1) immature granulocytes/bands (another type of immature white blood cells found in infections); 2) variant lymphocytes (an "atypical" circulating lymphocyte  
20 found in viral infections such as infectious mononucleosis); and 3) red blood cell (RBC) and platelet abnormalities such as RBC fragmentation, RBC agglutination, circulating nucleated RBCs (NRBCs - which are normally found only in the bone marrow)  
25 and platelet clumps (platelets are the cellular elements which help the blood to clot).

-25-

In a preferred aspect, the system is implemented on an IBM® compatible computer operating under Microsoft Windows™ 3.1 software, available from Microsoft Corporation of Redmond, Washington, to provide an appropriate user interface and interaction for a technician or physician. The programming of the modules 12a, 12b and 12c can be implemented in a conventional manner once knowledge of the invention is known from the description herein, for example, by the use of the commercially available Turbo Pascal™ for Windows™ software, and the programming of databases 14 can be implemented in a conventional manner with the commercially available Borland Paradox™ Engine 3.01, software, all of which are well known to those of ordinary skill in the art.

With respect to the database 14, the structure of each table in the database in the preferred implementation is shown in Appendices 1, 2, 3, and 4, shown at the end of this specification, for database components 1) patient demographics; 2) hemogram/peripheral blood data; 3) FCM data; and 4) bone marrow data; respectively.

Each module 12a, 12b,...12n includes a knowledge base 13 as defined previously. The knowledge base 13 of each of the decision support software modules 12a, 12b,...12n, for example, the

-26-

hematology module 12a, contains definitions of each of the patterns such that more than one pattern can coexist in the same case. For example, in acute leukemia, a thrombocytopenia pattern and/or normochromic, normocytic anemia pattern maybe present in addition to an abnormal mononuclear cell pattern. Since the goal is to arrive at a single dominant pattern, hierarchical classification of the patterns is arranged based upon the principles that, if WBC abnormalities are present, associated RBC and platelet abnormalities tend to be secondary to the WBC abnormalities, and the dominant pattern is the one that suggests a more specific, i.e., shorter, list of differential diagnoses. The latter rule is inferred from a commonly accepted heuristic criterion of diagnostic reasoning known as "Occam's razor" which states that the least complex of competing hypotheses is more likely to depict reality than the more complex hypotheses.

Thus, in implementation the hematology analyzer 18 values are displayed on the display 32, for example, a flat-panel display or a "cathode ray tube" (CRT) display, and referred to by the user of the device, e.g., a technician generally designated as box 22 of FIG. 1. Illustrative of what is displayed and referred to by the technician 22 on the display 32 are FIGS. 5 and 8. FIG. 5

-27-

illustrates the first screen displayed by the hematology module 12a which shows the hemogram parameters received from the hematology analyzer, including a suspect flag for NRBCs. FIG. 8 illustrates the hemogram findings received from the hematology analyzer including a suspect flag for blasts which are immature white blood cells that are only found circulating in the blood in large numbers when the diagnosis is acute leukemia or a related condition.

The hematology module 12a then automatically calculates a preliminary pattern based on the hemogram findings when the last analyzer parameter is received and the age and sex of the patient are known. The preliminary pattern is then displayed in a scrollable window of the display 32 and the hematology module 12a recommends a specific approach to the review of the peripheral blood smear based on the patterns that were present in the conflict set after the first pass through the data. This is illustrated by FIGS. 6 and 9.

FIG. 6 illustrates the second screen generated by hematology module 12a in which the system has indicated the NRBC suspect flag by automatically setting the check box "NRBCs." In this illustration the system has made a preliminary interpretation (based upon hemogram findings alone)

-28-

which is then displayed in a scrollable field shown in the upper half of FIG. 6.

Because of the NRBC suspect flag, the system has automatically set the "accept automated differential?" radio button to "No" and has indicated a morphologic review of white blood cells is necessary by setting the WBC check box. In FIG. 9 the preliminary indication generated by hematology module 12a includes a recommendation for reviewing RBC and WBC morphology and a manual WBC differential. The technologist 22 has agreed with the system recommendations after viewing the peripheral blood smear at scanning magnification.

As shown in FIG. 1, the technologist or user as illustrated by block 22 has complete control and has access to an electronic textbook 16 which allows access to didactic information as illustrated by FIGS. 20 and 21. The electronic textbook can be implemented, for example, in the commercially available electronic document design and publishing software product called Guide™, which is available from Infoaccess, Inc. of Bellevue, Washington. In FIG. 20 there is shown the "Table of Contents" window of the textbook. FIG. 21 illustrates a screen from within the textbook in the chapter on the "Leukoerythroblastic" pattern. The text is displayed in "Hypertext" format. The boldface text

-29-

with the arrow symbol is displayed in blue on the computer terminal indicating that the heading can be expanded to show more information. In FIG. 21, the user or technician 22 has clicked on the underlined  
5 text "except in the neonatal period" which brings up a comment box with further information about NRBCs in the newborn infants.

Thus, for example, in a case where the hematology analyzer 18 in FIGS. 1 and 2 reports a  
10 false, positive flag for NRBCs on the basis of the NRBC flag, for example, as illustrated in FIG. 5, the hematology module 12a would calculate and recommend a manual WBC differential for two reasons:  
1) the WBC count may need to be adjusted if a  
15 significant number of NRBCs were counted as leukocytes; and 2) the presence of circulating NRBCs may be associated with immature myeloid precursors in the peripheral blood. If the technologist 22, after scanning the peripheral blood smear, concludes  
20 that there are no circulating NRBCs, the technologist 22 notes that fact by un-highlighting the NRBC check box on the system. If the NRBC flag is the only feature which is causing the hematology module 12a to recommend a manual WBC differential,  
25 the system responds by withdrawing its recommendation for a manual differential. This is

-30-

illustrated by FIG. 7 which is the same screen as FIG. 6 showing changes after interaction with the technologist 22.

The technologist 22, after scanning the peripheral blood smear, has performed two functions as follows: 1) no NRBCs were visible on the peripheral blood smear (i.e., the analyzer suspect flag for NRBCs was a false positive). Therefore, the technologist 22 has unchecked the NRBC check box by clicking on the box. The system responded by automatically resetting the "Accept Automated Differential?" radio button to "Yes" since the only hemogram finding which supported a manual WBC differential was the NRBC suspect flag; and 2) the technologist 22 has also unchecked the WBC check box indicating that WBC morphology was normal at scanning magnification.

Likewise, the system 10 may recommend additional types of analysis such as, for example, a flow cytometry test from a flow cytometer 20 in which a flow cytometry module 12b comes into play with its own set of patterns (Table II) stored within the corresponding knowledge base 13 to be combined in a hierarchical manner to compliment the pattern operated on by the hematology module 12a. As may be appreciated, it is also possible to add additional tests as illustrated by FIG. 12, such as



-31-

bone marrow testing which results are processed by a bone marrow module 12c with its own patterns (Table III), the results of the test being reviewed by a physician, generally shown as box 24 in FIG. 1, as well as other tests as generally shown by the box labeled "XXX" and identified as 12n in FIG. 1. Such tests can include, for example, lymph node testing which results, as in the case with bone marrow testing, are all manually input through the keyboard 34 of FIG. 2.

In FIG. 12 there is illustrated the final interpretation and recommendation screen of the hematology blood module 12a, known as "Professor Petrushka." The highest pattern matched in the hierarchy "abnormal mononuclear cell pattern" is displayed in a scrollable field along with: 1) the specific findings in the peripheral blood; 2) the differential diagnosis (which, in this case, consists of only one diagnosis); and 3) a statement as to the specific diagnosis "acute myeloid leukemia" suggested by the finding "auer rods present." The recommendations for follow-up tests in this case are shown as including a bone marrow aspirate/biopsy, cytochemical stains, and immunophenotyping by flow cytometry.

The interaction between the modules 12a, 12b and 12c results in a multi-parameter report 26 which

-32-

results in a recommended diagnosis. In the case of the technologist 22 as well as the physician 24, it is also possible to access the electronic textbook 16 which explains the different conditions or patterns indicated by the hematology module 12a, flow cytometry module 12b and/or bone marrow module 12c to allow greater understanding of the diagnostic process as illustrated by FIGS. 20 and 21. The electronic textbook, as in the case with the database, would also be stored in DASD 46 with the knowledge base 13 and the modules 12a, 12b, 12c,...12n.

The electronic textbook 16 can be constructed from information obtained from the texts: "Wintrobe Clinical Hematology," 9th Edition, G.Richard Lee, Thomas C. Bithell, John Foerster, John W. Athens and John N. Lukens (editors), Lea & Febiger, Philadelphia, 1993; "Hematology of Infancy and Childhood," 4th Edition, David G. Nathan and Frank A. Oski (editors), W.B. Saunders, Philadelphia, 1992; and "Blood Textbook of Hematology," James H. Jandl, Little, Brown and Company, Boston 1987. As previously noted, the substance of these texts can be implemented into the electronics textbook 16 in a conventional manner by those of ordinary skill in the art, with, for example, the electronic document design and publishing software available under the

-33-

name Guide™ from InfoAccess, Inc., of Bellevue, Washington.

Having generally described the system and environment of the invention, FIG. 22 more clearly  
5 illustrates the method of the invention in a flow chart 50.

In its present implementation, the system in accordance of the invention consists of three knowledge based modules, i.e., a hematology module  
10 12a, a flow cytometry module 12b and a bone marrow module 12c, along with associated databases 14, an electronic textbook 16 with digitized images of, for example, abnormal blood cells, and in addition, educational materials available from Coulter  
15 Corporation of Miami, Florida under the name "Coulter VCS Technology Case Studies," volumes 1 and 2, 1992, showing case studies which are entered in the same format as the information from the previously mentioned textbooks. The modules 12a,  
20 12b and 12c serve to interpret data from a hematology blood analyzer such as a Coulter® STKS analyzer, a flow cytometry analyzer for interpretation of flow cytometry data, i.e., immunophenotyping and DNA content analysis, such as  
25 a Coulter® EPICS series flow cytometer, and bone marrow analysis of bone marrow morphology.

-34-

As shown in FIG. 22, there are fourteen steps implemented by the system modules 12a, 12b and 12c. Only the hematology module 12a follows all fourteen steps with the calculation of a preliminary pattern.

5 The flow cytometry module 12b goes directly from step 56 to step 68 and displays the final interpretation without user interaction. The bone marrow module 12c obtains the necessary peripheral blood data from the database 14 in step 54 and  
10 displays it along with suggestions for additional tests, such as, flow cytometry, i.e., FCM, special stains, cytogenetics, etc., based on the peripheral blood, i.e., PB, pattern. Then the bone marrow module goes directly from step 54 to step 66,  
15 performing the functions of step 56, i.e., data abstraction as illustrated hereafter, after all the data is collected, i.e., between steps 66 and 68. All three modules 12a, 12b and 12c, have incorporated therein the calculations of steps 56,  
20 58 and 60, except that all of the data is taken into consideration at that point. The term "data abstraction" refers to the process of comparing the raw data to the laboratory's normal ranges and the

-35-

medically important values that have been obtained from, for example, the texts employed to build the electronic textbook 16, resulting in a series of input values to step 58.

5           As previously discussed, knowledge representations in all three systems are based on defined patterns which serve as intermediate hypotheses in the diagnostic reasoning process. Within each pattern, there are additional "rules"

10   obtained from standard textbooks of hematology as discussed previously with reference to establishing the electronic textbook 16, which are followed to generate the text to the final interpretations and recommendations. This type of rule, an "if...then"

15   statement, e.g., "if observation 1 [and/or observation 2]... then conclusion," is well known to those of ordinary skill in the art. Specific examples of how these "rules" can be implemented to generate a text report are illustrated in Appendices

20   5, 6 and 7, which contain example code, and are shown at the end of this specification for the hematology module 12a, flow cytometry module 12b and bone marrow module 12c, respectively.

          It can be appreciated that an important

25   feature which distinguishes the present invention from other systems, is the method of knowledge representation implemented by the system. In

-36-

accordance with the invention there is first established a list of predetermined patterns, i.e., the names of each of the patterns. Then, the knowledge about the patterns is fixed by answering  
5 the following three questions about each pattern:

1) what minimum set of input data must be present to establish a positive match with the pattern under consideration?

2) is there any single piece of input data,  
10 or combination of input data, which, when present, rules out, i.e., excludes, that pattern from further consideration? and

3) what other patterns on the predetermined list are lower on the hierarchy than the pattern  
15 being considered, i.e., what other patterns can be "swallowed" (hereafter referred to as "swallowing" or "swallowing mechanism" which refers to the arrangement of patterns on the basis of which patterns predominate others) by the pattern under  
20 consideration?

These questions are answered in the practice of the invention, as shown from the tables, and the description of the program steps of FIG. 22 as discussed hereafter.

25 A description of the flow diagram 50 of FIG. 22 follows.

-37-

At step 52, a hematology workstation, for example, consisting of a IBM®-compatible computer with associated peripherals as shown in FIG. 3 and the software including modules 12a, 12b, 12c,...12n is interfaced to a hematology analyzer 18 and a flow cytometer 20, for example, such as the commercially available Coulter® STKS hematology analyzer and commercially available Coulter® EPICS Series flow cytometer, both available from Coulter Corporation of Miami, Florida. An interface 36 directly connects the hematology analyzer 18 to the computer 30. As data is transmitted from the hematology analyzer 18, alternatively from a flow cytometer 20, or input by hand through a keyboard 34 as shown in FIG. 13 which illustrates the screen dialog box for manually inputting the results of the flow cytometer 20 into the flow cytometry module (FCM) 12b, it is stored in a set of relational databases. As previously noted, such databases can be implemented through programming created through the use of the commercially available Borland Paradox™ Engine software program (see Appendices 1-4 showing example code). In addition to the patient's medical record number, age, sex and date of birth, the demographics database has fields for clinical history and other pertinent laboratory tests related to hematologic disease. The clinical information is entered into

-38-

the system manually and is stored in the demographics database (see Appendix 1).

The hemogram/peripheral blood (PB) database stores all of the numerical parameters normally transmitted from a hematology analyzer 18, as well as the suspect flags for morphologic abnormalities in the RBC, WBC, and platelet populations transmitted by the hematology analyzer (appendix 2). Wherever, in Appendices 1 through 4, the field type is referred to as "alphanumeric" (coded bitmap), this means that the data is coded in a string of characters.

For example, in Appendix 2, item 29, "suspect flags" is listed as field type alphanumeric (coded bitmap). There are 8 suspect flags transmitted by the hematology analyzer which are stored in this field. These flags are: 1) blasts; 2) immature granulocytes/bands; 3) variant lymphs; 4) NRBCs; 5) dimorphic RBC; 6) micro RBCs/RBC fragments; 7) RBC agglutination; and 8) platelet clumps (see FIG. 5). Field 29 of the hemogram database, "suspect flags," is a string of 8 characters, conventionally numbered characters 0 through character 7 using the nomenclature of so-called "null-terminated strings" known to individuals of ordinary skill in the field of computer programming. The result of the suspect flag "blasts," i.e., present or absent, is stored in



-39-

character 0 of the string as a "0" (absent) or "1" (present); the result of the flag "immature granulocytes/bands," i.e., "absent," "flag type 1," or "flag type 2," is stored in character 1 of the string as a "0" (absent), "1" (Flag type 1) or "2" (flag type 2). The remaining flags can only be present or absent, and are stored in the subsequent characters of the string exactly as for the "blasts" flag. Therefore, if the hematology analyzer transmits only a "variant lymph" flag and all of the other flags are negative, the character string "00100000" is stored in field 29 "suspect flags."

In step 54, the operator chooses a case, i.e., a patient/sample, and the system displays data retrieved from the database 14. The hematology module 12a displays the hemogram findings, the flow cytometry module 12b displays all of the flow cytometry data.

FIG. 14 shows the screen from the flow cytometry module 12b displaying the flow cytometry data on a lymph node biopsy in a patient with "chronic lymphocytic leukemia" (CLL). The data (displayed on the left) were manually input by the user as fluorescence intensity. FIG. 15 shows the screen from the flow cytometry module 12c showing the interpretation of flow cytometry testing performed on a peripheral blood specimen from a

-40-

patient with "acute lymphoblastic leukemia" (ALL).  
The flow cytometry data was input by the user as a  
percentage of positive cells.

The bone marrow module 12c looks up the  
5 patient's blood findings in the hemogram/PB database  
and displays them along with recommendations for the  
bone marrow which are triggered by the PB pattern.  
FIG. 16 shows the first screen in the bone marrow  
module 12c when the user entered the specimen  
10 identification number into a dialog box (not shown),  
the system 10 responded by looking up the peripheral  
blood information on this patient saved in the  
database 14. The peripheral blood findings are  
displayed along with recommendations as to how to  
15 approach the bone marrow examination. In this case  
there are reminders for the additional tests which  
need to be performed on the bone marrow.

The type of program to be implemented for  
retrieving such information from a database and  
20 displaying such information on the screen will be  
readily apparent to those of ordinary skill in the  
art of computer programming based on this  
description.

In step 56 all of the modules 12a-12c go  
25 through a process of data abstraction.  
Specifically, as previously noted, "data  
abstraction" refers to the process of comparing the

-41-

raw data, i.e. the instrument parameters, to the normal ranges for these parameters set by the laboratory, and to the medically important ranges for these parameters which are obtained from the standard hematology textbooks such as have been identified with reference to the electronic textbook 16. The concept of "data abstraction" and "heuristic classification" are well known, as is apparent from the text "Heuristic Classification," Clancey, W.J., "Artificial Intelligence," 1985; 27:289-350. The comparison between the instrument parameters and the normal and medically important ranges results in a set of "input values" to the computer program, i.e., the inference engine 15 designed as modules 12a, 12b, 12c,...12n. The input parameters for the peripheral blood module 12a, are listed in Table IV below.

TABLE IV

## PERIPHERAL BLOOD INPUT PARAMETERS

20	1. Age      1a=0-1 (infant) 1b=2-13 (child) 1c=14-35 (young adult) 1d=36-49 (adult) 1e=50-69 (middle age) 1f=greater than 69 (elderly)
	2. Sex      2a=male 2b=female
25	Actual numeric values for parameters 3-13 set by hospital/laboratory normal ranges and medically important ranges

-42-

- |    |   |  |
|----|---|--|
|    | 3. Hgb  | 3a=severe anemia 3b=moderate anemia<br>3c=mild anemia 3d=normal 3e=mild<br>polycythemia 3f=severe polycythemia |
|    | 4. MCV  | 4a=microcytosis 4b=normocytic 4c=macrocytosis  |
| 5  | 5. MCH  | 5a=hypochromic 5b=normochromic   |
|    | 6. RDW  | a=normal 6b=mild anisopoikilocytosis<br>6c=marked anisopoikilocytosis  |
|    | 7. PLT  | 7a=thrombocytopenia 7b=normal<br>7c=thrombocytosis   |
| 10 | 8. WBC  | 8a=severe leukopenia 8b=mild leukopenia<br>8c=normal 8d=mild leukocytosis 8e=severe<br>leukocytosis            |
|    | 9. Neut   | 9a=granulocytopenia 9b=normal<br>9c=granulocytosis   |
| 15 | 10. Lymph                                       | 10a=lymphocytopenia 10b=normal 10c=mild<br>lyphocytosis 10d=severe lymphocytosis                               |
|    | 11. Mono  | 11a=normal 11b=absolute monocytosis  |
|    | 12. Eos   | 12a=normal 12b=eosinophilia  |
|    | 13. Baso  | 13a=normal 13b=basophilia  |
| 20 | 14. WBC Population Flag                         | a=normal b=abnormal  |
|    | 15. RBC Population Flag                         | a=normal b=abnormal  |
|    | 16. PLT Population Flag                         | a=normal b=abnormal  |
|    | 17. Suspect Blasts Flag                         | a=present b=absent   |
| 25 | 18. Suspect Immature<br>Granulocytes/Bands Flag | a=present b=absent   |
|    | 19. Variant Lymphs Flag                         | a=present b=absent   |
|    | 20. NRBCs Flag                                  | a=present b=absent   |
|    | 21. MicroRBC/Fragments Flag                     | a=present b=absent   |
|    | 22.* Platelet Clumps Flag                       | a=present b=absent   |

-43-

23. Dimorphic RBC Flag            a=present b=absent  
 24.\* RBC Agglutination           a=present b=absent  
 25.\* MCHC                         Normal vs Increased

\*=Used only for setting a check box on the screen or  
 displaying text. Not considered by the program.

5

For the hematology module 12a there are up to 25  
 parameters, of which 22 are used by the program in  
 step 58 (see Table IV). For example, for parameter  
 #3 (Hgb), the results of "data abstraction" are as  
 follows, depending on the laboratory's normal range:  
 3a = severe anemia, 3b = moderate anemia, 3c = mild  
 anemia, and 3d = normal Hgb, 3e = mild polycythemia,  
 3f = severe polycythemia. For the flow cytometry  
 module 12b there are thirty-seven immunologic  
 markers and several DNA-analysis parameters such as  
 cell size, DNA index and percentage of cells in S-  
 phase which get abstracted and translated into  
 "input values" to the program. These parameters are  
 listed in Table V below.

10

15

20

TABLE V

## FLOW CYTOMETRY INPUT PARAMETERS

25

- |                  |          |
|------------------|----------|
| 1. Age           | 24. CD45 |
| 2. Sex           | 25. CD56 |
| 3. Specimen Type | 26. CD57 |
| 4. CD1a          | 27. CD61 |

-44-

5	5.	CD2	28.	CD71
	6.	CD3	29.	Kappa
	7.	CD4	30.	Lambda
	8.	CD5	31.	Ia (HLA-DR)
	9.	CD7	32.	Glycophorin A
10	10.	CD8	33.	Cytoplasmic Kappa
	11.	CD10	34.	Cytoplasmic Lambda
	12.	CD11c	35.	Cytoplasmic Mu
	13.	CD13	36.	Cytoplasmic CD3
	14.	MY4 (CD14)	37.	Cytoplasmic CD22
15	15.	Mo2 (CD14)	38.	Tdt
	16.	CD16	39.	Myeloperoxidase (by cytochemistry)
	17.	CD19	40.	Myeloperoxidase (by monoclonal antibody)
	18.	CD20	41.	Cell Size [Small, Medium, Large, Variable, Bimodal]
	19.	CD22	42.	DNA-Index
20	20.	CD25	43.	% S-phase
	21.	CD33	44.	% Abnormal Cells
	22.	CD34	45.	% Blasts or Promyelocytes Present
	23.	CD41	46.	Peripheral blood and clinical findings [from database 14]

For item 3, Specimen Types are: PB=Peripheral Blood, BM=Bone Marrow, LN=Lymph Node, SP=Spleen, FL=Fluid, EN=Extranodal

30 Items 4-40 are expressed either: (1) in percent positive (>30=positive); or, (2) by fluorescence intensity [-/+ ++/+++] (except item 39)

During Data Abstraction (Step 56) Items 4-40 are designated: POSITIVE (%>30 or + or ++ or +++), NEGATIVE (%<=30 or -) or NOT DONE.

35 For the bone marrow module 12c, the results of the user's input for the bone marrow differential and the qualitative/quantitative features from the aspirate and biopsy (cellularity, presence of abnormal cells, results of special stains, etc.)

-45-

make up the "input values." A complete listing of the bone marrow input parameters is found in Table VI below.

TABLE VI

5	BONE MARROW INPUT PARAMETERS
	1. Age
	2. Sex
	3. Specimen(s) [Aspirate/Biopsy/Touch Prep]
	4. Peripheral Blood Findings [from database 14]
10	5. Spicules [Present/Absent]
	6. Megakaryocytes [Normal/Decreased/Increased]
	7. Bone Marrow Cellularity [%]
	8. Abnormal Mononuclear Cells/Foreign Cells [Present/Absent/Type]
15	9. Bone Marrow Differential
	10. Erythroid Karyorrhexis [Present/Absent]
	11. Hypogranulation [Present/Absent]
	12. Abnormal Eosinophils [Present/Absent]
	13. Nuclear/Cytoplasmic Asynchrony [Absent/Myeloid/Erythroid/Myeloid and Erythroid]
20	14. Giant Precursors [Absent/Myeloid/Erythroid/Myeloid and Erythroid]
	15. Cytoplasmic Vacuoles [Absent/Myeloid/Erythroid/Myeloid and Erythroid]
25	16. Organisms
	17. Bony Trabeculae [Normal/Thick/Thin]
	18. Collagen Fibrosis [Present/Absent]
	19. Pattern of Lymphocyte/Plasma Cell Infiltration [Diffuse/Nodular/Paratrabecular]
30	20. Granulomas [Present/Absent]
	21. Iron Storage [Normal/Decreased/Increased]
	22. Iron Incorporation [Normal/Decreased/Increased]
	23. Ring Sideroblasts [%]
	24. Peroxidase Stain [Positive/Negative]
35	25. Non-Specific Esterase Stain [% positive]
	26. PAS Stain [Positive/Negative]
	27. TRAP Stain [Positive/Negative]
	28. Reticulin Stain [Negative/1+/2+/3+]
	29. Flow Cytometry Findings [from database 14]

-46-

For the hematology module 12a, the program in step 58 works as follows: The first 17 peripheral blood patterns, i.e., all of the patterns except for the "non-specific" pattern, are assembled into a set. At this point, at the start of the program, the set "patterns matched" indicated in computer programming by elements between square brackets contains 17 values as follows [1 2 3 4 5 6 7 8 9 10 11 12 13 14 15 16 17]. The input values from step 56 are then compared to the pattern definitions for each of the 17 patterns beginning with pattern 1 and continuing sequentially until pattern 17. The pattern definitions for the hematology patterns are shown in Table VII as follows.

15

TABLE VII

## PERIPHERAL BLOOD PATTERN DEFINITIONS

	<u>Pattern</u>	<u>Definition</u>
	1. Normal	No abnormality, No STKS suspect flags
20	2. Reactive Neutrophilia	Increased WBC Count (<50,000) and increased neutrophils
	3. Leukoerythroblastic	NRBC flag, Immature gran/band flag, Increased RDW
	4. Abnormal mononuclear cell	Blast flag or variant lymphocyte flag
25	5. Erythrocytosis	Increased Hgb



-47-

	6. Thrombocytosis	Increased platelet count
	7. Absolute lymphocytosis	Increased lymphocyte count and lymphocytes $\geq 35\%$
	8. Thrombocytopenia	Decreased platelet count
5	9. Eosinophila	Increased eosinophil count
	10. Monocytosis	Increased monocyte count and monocytes $\geq 15\%$
	11. Bi-pancytopenia	Decreased (two or more) of: Hgb, WBC count, Platelet count
10	12. Leukopenia	Decreased WBC count or Decreased Neutrophil count
	13. Hypochromic and/or	Decreased MCH or Decreased MCV microcytic
	14. Macrocytosis	Increased MCV
15	15. Normochrmic, normocytic anemia	Decreased Hgb, normal MCH, normal MCV
	16. Anemia with anisocytosis	Decreased Hgb, Increased RDW
	17. Leukemoid	Increased WBC count ( $>50,000$ ), Immature gran/band flag, Neut $\geq 65\%$
20	18. Non-specific	Anything which does not meet definitions 1-17

The pattern definitions are expressed within the program as simple "if...then" statements well known to those of ordinary skill in the art. The system in accordance with the invention eliminates those patterns whose definitions are not met. Therefore, the "if...then" statements are of the form, "if the pattern definition for pattern x is

-48-

not met, then eliminate pattern x from the set patterns matched." For example, it can be seen from Table IV above that the input value "7c" means thrombocytosis which is the medical term for an increased platelet count. It can also be seen from previously set forth Table VII above that the only parameter in the definition for pattern 6 (thrombocytosis) is an increased platelet count, i.e., thrombocytosis. In English sentence structure, the "if...then" statement for pattern 6 is therefore, "if the input value 7 is not 7c then remove pattern 6 from list of patterns matched." As mentioned above, each pattern from 1 to 17 is evaluated in exactly this manner according to its own definition.

The FCM module 12b and the bone marrow module 12c perform pattern matching in an exactly analogous manner but they do so in step 68 (final pattern) only. There is no preliminary pattern for the FCM and bone marrow modules 12b and 12c. The pattern definitions for the FCM module are shown in Table VIII below and the possible input parameters to step 68 for the FCM module 12b are listed in Table V above.

-49-

## VIII

## FLOW CYTOMETRY PATTERN DEFINITIONS

	<u>Pattern</u>	<u>Definition</u>
	1. Stem Cell	Positive for Tdt or CD34
5	2. Myeloid	Positive for Myeloperoxidase or CD13 or CD33
	3. Erythroid	Positive for Glycophorin A
10	4. Megakaryocytic	Positive for CD41 or CD61
	5. Immature B-cell	Positive for (Tdt or CD34) and (CD19, CD22,cCD22 or cmu)
	6. Immature T-cell	Positive for (Tdt or CD34) and (CD7 or cCD3)
15	7. Mature B-cell	Positive for Kappa or Lambda or (Tdt negative with CD19 or CD20 or CD22)
	8. Mature T-cell	CD3 or cCD3 or (CD2 and (CD2 or CD7 or CD8)) or (CD7 and CD4) or (CD7 and CD8)
20	9. Plasma cell	Positive for cKappa or cLambda
	10. Non-specific	(1) CD14 alone or (2) CD19 positive and Tdt, kappa, lambda ND or (3) CD3 or CD7 or cCD3 positive and Tdt, kappa, lambda ND or (4) Does not meet any definition 1-9.
25		
	<b>Abbreviations:</b> Tdt=Terminal deoxynucleotidal transferase; cCD3=cytoplasmic CD3; cCD22=cytoplasmic CD22; cKappa=cytoplasmic Kappa; cLambda=cytoplasmic Lambda ND=Not Done	
30		

-50-

For the FCM module 12b, during the data abstraction step, each of the immunophenotype input values are considered "Positive," "Negative" or "Not Done." Similarly, the pattern definitions for the bone marrow module are listed in Table IX and the input values are listed in Table VI above.

**TABLE IX**  
**BONE MARROW PATTERN DEFINITIONS**

	<u>Pattern</u>	<u>Definition</u>
10	1. Normocellular	Cellularity >30% and <70%, Blasts <5%, Lymphocytes <20%, Plasma cell <10%, Blasts <Promyelocytes + Myelocytes + Metamyelocytes, Promyelocytes + Myelocytes + Metamyelocytes <Neutrophils
15		No myeloid or erythroid predominance No qualitative abnormalities
	2. Hypocellular	Cellularity <30%
20	3. Erythroid predominance/ hyperplasia	Sum of erythroid precursors >= Sum of myeloid precursors
	4. Hypercellular, predominance/ hyperplasia	Cellularity >=70%, Myeloid > Erythroid, Sum of Myeloid >50%
25	5. Megakaryocytic hyperplasia	Megakaryocytes increased
	6. Lymphocytosis/ Plasmacytosis	Cellularity >30%, Lymphocytes >20% or Plasma cells >10%
30	7. Mononuclear Cell Infiltration	Blasts (non-erythroid) >=30% or Promyelocytes >=30% or Other abnormal mononuclear cells present

-51-

- |   |                              |  |
|---|------------------------------|--|
| 5 | 8. Foreign Cell Infiltration | Foreign cells present                                      |
|   | 9. Granulomatous             | Granulomas present   |
|   | 10. Bi-hyperplasia           | (Patterns 4 and 5 together) or (Patterns 3 and 5 together) |
|   | 11. Non-specific             | Anything which does not meet any definition 1-10.          |

In step 60, the patterns which remain after step 58 are then considered. Considering again the peripheral blood module 12a, if all of the specific patterns have been eliminated in step 58, i.e., in mathematical and computer terms, only the empty set, designated [ ], remains, then the non-specific pattern, i.e. pattern 18, is added to the set. In this case, the non-specific pattern becomes the predominant pattern by default. If there was one specific pattern left in the set after step 58, it becomes the predominant pattern. If there is more than one specific pattern left after step 58, the module 12a branches to the "swallowing mechanism," previously described. As noted, this is one of the most important steps in the process and is based on knowledge obtained from standard textbooks of hematology such as those in the electronic textbook 16 previously discussed. This process is based on the established hierarchy of the patterns. Specifically, the patterns remaining in the set of

-52-

"patterns matched" which are swallowed by other patterns also remaining in the set are listed in Table X below.

TABLE X

5	PERIPHERAL BLOOD SWALLOWING	
	<u>Pattern</u>	<u>Swallows Pattern Numbers</u>
	1. Normal	NONE
	2. Reactive Neutrophilia	6 (except Plt > 1,000,000), 9, 10, 13, 14, 15, 16
10	3. Leukoerythroblastic	2, 5, 6, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17
	4. Abnormal mononuclear cell	2, 3, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17
	5. Erythrocytosis	2, 6, 9, 10, 13, 14
15	6. Thrombocytosis	9, 10 (if Plt > 1,000,000)
	7. Absolute Lymphocytosis	if (Age > 45 or Severe) then 2, 3, 5, 6, 8, 9, 10, 11, 12, 13, 14, 15, 16 else 9, 10
	8. Thrombocytopenia	2, 5, 9, 10, 13, 14, 15, 16
20	9. Eosinophilia	15, 16
	10. Monocytosis	6 (if Plt 500,000-999,000), 9, 12, 13, 14, 15, 16
	11. Bi-pancytopenia	2, 5, 6, 10, 13, 14
25	12. Leukopenia	5, 6 (except Plt > 1,000,000), 8, 9, 13, 14, 15, 16, 17
	13. Hypochromic and/or microcytic	6 (except Plt > 1,000,000), 9, 14, 15, 16
	14. Macrocytosis	6, 9, 13, 15, 16

-53-

- |                                     |  |
|-------------------------------------|--|
| 15. Normochromic, normocytic anemia | 6 (except Plt>1,000,000)                 |
| 16. Anemia with anisocytosis        | 6 (except Plt>1,000,000), 15             |
| 17. Leukemoid                       | 2, 5, 6, 8, 9, 10, 11, 13, 14, 15,<br>16 |

5 For example, if the remaining set contains patterns  
4, 8 and 15, then the swallowing mechanism in step  
60 eliminates both patterns 8 and 15 from the set  
because, as can be seen in Table X, pattern 4  
"swallows" all of the other patterns, i.e., it is at  
10 the top of the hierarchy. In this example, after  
going through the swallowing mechanism in step 60,  
pattern 4 becomes the predominant pattern.

The specific rules for the swallowing  
mechanism in step 60 for the FCM module 12b and the  
15 bone marrow module 12c are given in Tables XI and  
XII respectively, as set forth below.

TABLE XI

## FLOW CYTOMETRY SWALLOWING

20	<u>Pattern</u>	<u>Swallows - "Patterns Not Excluded With Numbers"</u>
	1. Stem Cell	7, 8, 9
	2. Myeloid	1, 7, 8, 9
	3. Erythroid	1, 5, 6, 7, 8, 9
	4. Megakaryocytic	1, 7, 8, 9
25	5. Immature B-cell	1, 3, 4, 7, 8, 9
	6. Immature T-cell	1, 3, 4, 7, 8, 9
	7. Mature B-cell	1, 2, 3, 4, 5, 6, 8, 9

- 54 -

- |                  |                         |
|------------------|-------------------------|
| 8. Mature T-cell | 1, 2, 3, 4, 5, 6*, 7, 9 |
| 9. Plasma cell   | 1, 2, 3, 4, 5, 6, 7, 8  |
| 10. Non-specific | NONE                    |

\* if CD4+/CD8- or CD4-  
/CD8+

5

TABLE XII

## BONE MARROW SWALLOWING

	<u>Pattern</u>	<u>Swallows Pattern Numbers</u>
	1. Normocellular	NONE
10	2. Hypocellular	1, 3, 5
	3. Erythroid predominance/hyperplasia	1
	4. Hypercellular, myeloid predominance/ hyperplasia	1
	5. Megakaryocytic hyperplasia	1
15	6. Lymphocytosis/plasmacytosis	1, 3, 4
	7. Mononuclear cell infiltration	1, 2, 3, 4, 5, 6, 9, 10
	8. Foreign cell infiltration	1, 2, 3, 4, 5, 6, 9, 10
	9. Granulomatous	1, 10
	10. Bi-hyperplasia	1, 3, 4, 5

20           The flow cytometry module 12b does  
 "swallowing" (in step 68) a little bit differently  
 than the hematology module 12a and bone marrow  
 module 12c. In certain cases, two or more patterns  
 can occur simultaneously, e.g., a biphenotypic  
 25 leukemia can meet the definition of both the myeloid  
 pattern and the immature B-cell pattern. Instead of  
 swallowing patterns from the set "patterns matched,"  
 the flow cytometry module 12b keeps track of which



-55-

patterns have not been excluded because certain antibodies have not been tested, and compares that list (the set "patterns not excluded") with a list of the patterns that have already been positively  
5 matched. The positively matched patterns then "swallow" patterns from the "patterns not excluded" list depending on the hierarchy. In this way, if the test results are not conclusive, and the flow cytometry module 12b requests that additional  
10 antibodies be studied, the program only asks the user to test for additional possibilities that are logical, i.e., those that remain in the "not yet excluded" list after swallowing has taken place.

The bone marrow module 12c uses the same  
15 swallowing mechanism (i.e., swallowing of patterns from the set "patterns matched") as the peripheral blood module 12a.

In step 62, hematology module 12a compares the specific hemogram findings and the preliminary  
20 predominant pattern to a series of "rules" taken from the previously referenced standard textbooks of hematology for electronics textbook 16 which generate recommendations such as "no peripheral blood smear is necessary" or "scan the smear," "look  
25 at WBC morphology" and "perform a manual differential." The rules for generating this text are of a similar nature to the examples of

-56-

programming code in appendices 5, 6 and 7. Programming of this type using "if...then" statements is well known to those of ordinary skill in the art.

5           In step 64 hematology module 12a displays the pattern, the hemogram findings (equivalent to the commercially available Coulter® STKS hematology analyzer definitive flags) and recommendations calculated in steps 56, 58 and 60 in normal text to  
10 the user. The program also sets certain check boxes on the display 32 such as "RBC agglutination present" depending on the suspect flags indicated by, e.g., the previously referenced Coulter® Educational Series text (FIGS. 6 and 9).

15           In step 66, the user enters peripheral blood or blood marrow findings interactively with hematology module 12a or bone marrow module 12c, as respectively shown in FIGS. 7, 10, 11, 17 and 18.

FIG. 7 has been previously described. FIG.  
20 10 shows the RBC morphology screen where the technologist 22 has indicated "1 + polychromasia" and "1 + macrocytes" by clicking on appropriate radio buttons. FIG. 11 is the WBC screen. The  
25 technologist 22 has entered the manual differential into a dialog box (not shown) and the system has displayed that differential which included 86% blasts. The technologist 22 has indicated that

-57-

"auer rods" are present in the blasts by selecting "present" from the drop-down list box. Auer rods are abnormal granules that are only seen in blasts from patients with acute myeloid leukemia. FIG. 17 is the entry screen for bone marrow aspirate data. The physician 24 has indicated the findings "spicules present" and "megakaryocytes decreased" by clicking on the appropriate radio buttons. The bone marrow differential was entered into a dialog box (not shown) and is displayed by the system. FIG. 18 is the screen for inputting the results of special stains on the bone marrow. The physician 24 has input "peroxidase positive blasts" which indicates an acute myeloid leukemia and NSE (non-specific esterase) positive in 90% of the cells (this finding indicates a monocytic leukemia which is designated type M5a in the FAB, i.e., the French-American-British classification of acute myeloid leukemia). Thus, as shown, the user is free to follow the programs suggestions for additional testing (FIG. 9) or not follow them as they wish (FIG. 7). The module 12a or 12c keeps track of the user input and updates the algorithm "input values" automatically during this process. For example, if the STKS "suspect blast" flag was present, input value #17 is

-58-

set to "17a." If the WBC manual differential showed no blasts, the human input overrides the hematology analyzer input and parameter #17 is set to "17b."

In step 68 "data abstraction," comparison of  
5 the raw data with the established ranges for bone marrow module 12c and flow cytometry module 12b takes place. As mentioned above, in this step all of the modules 12a, 12b and 12c go through the pattern matching program including swallowing, if  
10 necessary, exactly as in hematology module 12a at steps 56, 58 and 60, except that all of the data is taken into consideration resulting in a final predominant pattern.

In step 70 all of the modules 12a, 12b and  
15 12c branch to a procedure which is a set of "diagnostic rules" for each given predominant pattern. The rules are used to generate the text of the final report, that is the module's final interpretation/recommendations. Examples of the  
20 rules to generate text, as previously noted, are given in Appendices 5, 6, and 7. Also, this type of programming uses the technique of "if...then" statements well known to those of ordinary skill in the art, and the rules themselves are taken from  
25 standard textbooks of hematology as previously referred to for the electronic textbook 16.

-59-

In step 72, in conjunction with step 70, and before assembling the final report, the system checks all of the appropriate databases for additional information. This is a very important step in the process. The databases are indexed on the basis of the patient's identification number. Hematology module 12a checks for certain items of clinical history in the demographics portion of database 14 and also for previous blood specimens on the same patient. The "diagnostic rules" within the final pattern, obtained from the referenced standard textbooks of hematology, govern how to use this information to modify the final report as shown illustratively in Appendices 5, 6 and 7.

Flow cytometry module 12b checks for history, blood and/or bone marrow findings, and previous FCM data. Notice that from FIG. 14, the box labelled "Additional Information" is empty. In this patient, the flow cytometry test was performed on the lymph node biopsy before the peripheral blood test was performed. For this reason, there was no peripheral blood information in the database 14 on this patient. In FIG. 15, the designation "L3" in the box "Additional Information" indicates that the system found peripheral blood data in the database 14 and that the technologist had found blasts with "Burkitt's cell morphology" (designated type L3 by

-60-

the "French-American-British" [FAB] classification of acute lymphoblastic leukemia) on the manual differential of the peripheral blood. The interpretation includes the pattern matched, the differential diagnosis, a correlation of the findings in the database 14 with flow cytometry data, and a recommendation to confirm these results with a bone marrow/aspirate biopsy. Bone marrow module 12c already knows about the current peripheral blood specimen which was retrieved in step 54. Bone marrow module 12c looks up the history, any FCM data (from flow cytometry module 12b's database) and any previous bone marrow report.

In step 74, before generating the final report, the peripheral blood module 12a goes through another procedure to check if any important information was suppressed by the "swallowing" mechanism. This is necessary because some of the findings may be the result of two pathophysiological processes instead of just one disease. If important information is found, this generates text that is called "drop-in" packages. The rules for the use of drop-in packages are specifically stated in Table XIII below.

-61-

TABLE XIII

## PERIPHERAL BLOOD DROP-IN PACKAGES

	If this pattern was "swallowed":	Then "drop-in" comments when final predominant pattern is:
5		
	2 (Reactive Neutrophilia)	8 (Thrombocytopenia) 11 (Bicytopenia but with increased WBC) 13 (Hypochromic and/or microcytic)
10		
	13 (Hypochromic and/or microcytic)	Any pattern except 3 (Leukoerythroblastic)
	14 (Macrocytosis)	2 (Reactive Neutrophilia) 8 (Thrombocytopenia) 12 (Leukopenia)
15		
	15 (Normochromic, normocytic anemia)	7 (Absolute Lymphocytosis) 8 (Thrombocytopenia)

For example, according to Table XIII, a macrocytosis drop-in package (pattern #14) will be added in a patient who has a reactive neutrophilia pattern (predominant pattern #2) because the macrocytosis (i.e., enlargement of the red blood cells) is likely to be due to a separate process (e.g., vitamin deficiency) than the one causing the increased number of white blood cells (e.g., an infection). That is, a patient with both macrocytosis and reactive neutrophilia might have both an infection

-62-

and a vitamin deficiency. The text for a drop-in package is a shortened version of the text generated when that pattern is the predominant pattern.

5 In step 76 the final interpretation and recommendations are displayed in a user editable field (a "hard copy" of this report can be generated if the laboratory wishes to do so, for example, at printer 39).

10 For hematology module 12a and flow cytometry module 12b, this "report" contains the patient identification, final predominant pattern (could be more than one in flow cytometry module 12b), the differential diagnosis, a suggested interpretation if the findings are specific enough, and  
15 recommendations for additional history or laboratory tests. This is shown in FIG. 14 for the flow cytometry module where the final interpretation of these data are displayed in the scrollable field. The interpretation includes: 1) the pattern  
20 matched, "mature B-cell pattern (monoclonal for lambda)"; 2) the differential diagnosis; 3) the findings which suggest the specific diagnosis of (CLL) or "small lymphocytic lymphoma"; and, 4) a recommendation to correlate these results with other  
25 available information. FIG. 15 shows correlation with additional information by the designation "L3" as previously discussed.



-63-

For bone marrow module 12c, the final report is a bone marrow report as a pathologist would dictate it, including peripheral blood information, a description of all the pertinent bone marrow findings, the results of the FCM immunophenotyping, and a final interpretation and comments. This is shown in FIG. 19 which illustrates the final interpretation screen for the bone marrow module 12c. The screen consists of a scrollable field in which is displayed the final bone marrow report. The patient identification information has been scrolled off at the top of the screen to show the remainder of the interpretation. The report contains: 1) the peripheral blood findings retrieved from the database 14; 2) the highest bone marrow pattern matched in the hierarchy, "mononuclear cell infiltration"; 3) the bone marrow morphologic findings entered by the physician; 4) a summary of the immunology interpretation (flow cytometry module 12c) retrieved from the database 14; and 5) the final interpretation with the FAB classification.

In step 78 the findings and interpretations are coded and saved to the appropriate portion of database 14. The information entered by the technologist 22 or physician 24 by interaction with the various Windows™ controls such as check boxes,

-64-

radio buttons and drop-down lists are coded in the bitmap scheme described above. There are alphanumeric fields in the database 14 corresponding to each of the screens which can accept user input in a module (e.g., FIG. 6) from the peripheral blood module 12a corresponds to the field #31 "scan of the peripheral blood smear" in the hemogram database (Appendix 2). These fields contains one character for each of the controls on the corresponding screen. The state of each of the controls when the user clicked on the "proceed" button (i.e., for check boxes and radio buttons whether or not it was checked and for drop-down lists the number of the item that was selected) is coded in the appropriate character for that field. As discussed previously throughout the specification, the following appendixes will aid in understanding the invention when viewed with references.

## APPENDIX 1: DEMOGRAPHICS DATABASE STRUCTURE

	<u>Field Name:</u>	<u>Field Type:</u>
	1. Patient Identification Number	Alphanumeric
	2. Patient Name	Alphanumeric
5	3. Sex	Alphanumeric
	4. Date of Birth	Paradox Date (long integer)
	5. Race	Alphanumeric
	6. Inherited Diseases	Alphanumeric (coded bitmap)
	7. Social History	Alphanumeric (coded bitmap)
10	8. Previous Surgery	Alphanumeric (coded bitmap)
	9. Drug History	Alphanumeric (coded bitmap)
	10. Date of Drug History	Paradox Date (long integer)
	11. Cancer History	Alphanumeric (coded bitmap)
	12. Autoimmune History	Alphanumeric (coded bitmap)
15	13. Endocrine History	Alphanumeric (coded bitmap)
	14. Impaired Immunity	Alphanumeric (coded bitmap)
	15. Lab Tests	Alphanumeric (coded bitmap)

-66-

## APPENDIX 2: HEMOGRAM DATABASE STRUCTURE

	<u>Field Name:</u>	<u>Field Type:</u>
	1. Date of Specimen	Paradox Date (long integer)
	2. Time of Specimen	Alphanumeric
5	3. Patient Identification Number	Alphanumeric
	4. Sequence Number	Alphanumeric
	5. Cassette/Position	Alphanumeric
	6. Patient Location	Alphanumeric
	7. Physician	Alphanumeric
10	8. User Field #1	Alphanumeric
	9. User Field #2	Alphanumeric
	10. User Field #3	Alphanumeric
	11. Comments	Alphanumeric
	12. RBC Count	Numeric
15	13. Hemoglobin	Numeric
	14. Hematocrit	Numeric
	15. MCV	Numeric
	16. MCH	Numeric
	17. MCHC	Numeric
20	18. RDW	Numeric
	19. Platelet Count	Numeric
	20. WBC Count	Numeric
	21. Neutrophils [%]	Numeric
	22. Lymphocytes [%]	Numeric
25	23. Monocytes [%]	Numeric
	24. Eosinophils [%]	Numeric
	25. Basophils [%]	Numeric
	26. WBC Population Flag	Alphanumeric
	27. RBC Population Flag	Alphanumeric
30	28. Platelet Population Flag	Alphanumeric
	29. Suspect Flags	Alphanumeric (coded bitmap)
	30. Preliminary Pattern	Numeric
	31. Scan of Peripheral Blood Smear	Alphanumeric (coded bitmap)
	32. RBC Morphology	Alphanumeric (coded bitmap)
35	33. Organisms	Alphanumeric
	34. Manual Differential -- Neutrophils	Numeric
	35. Manual Differential -- Bands	Numeric
	36. Manual Differential -- Lymphocytes	Numeric
	37. Manual Differential -- Monocytes	Numeric
40	38. Manual Differential -- Eosinophils	Numeric
	39. Manual Differential -- Basophils	Numeric
	40. Manual Differential -- NRBCs	Numeric
	41. Manual Differential -- Blasts	Numeric
	42. Manual Differential -- Promyelocytes	Numeric
45	43. Manual Differential -- Myelocytes	Numeric
	44. Manual Differential -- Metamyelocytes	Numeric
	45. Other Mononuclear Cells (OMNC)	Numeric

-67-

5	46. Type of OMNC	Alphanumeric
	47. Manual Platelet Count	Numeric
	48. WBC Morphology	Alphanumeric (coded bitmap)
	49. WBC Inclusions	Alphanumeric
	50. Final Pattern	Numeric
	51. Reticulocyte Count	Numeric

## APPENDIX 3: FLOW CYTOMETRY DATABASE STRUCTURE

	<u>Field Name:</u>	<u>Field Type:</u>
	1. Specimen Identification Number	Alphanumeric
	2. Patient Identification Number	Alphanumeric
5	3. Clinical Information	Alphanumeric
	4. Specimen Type	Alphanumeric
	5. Specimen Date	Paradox Date (long integer)
	6. CD1a	Alphanumeric
	7. CD2	Alphanumeric
10	8. CD3	Alphanumeric
	9. CD4	Alphanumeric
	10. CD5	Alphanumeric
	11. CD7	Alphanumeric
	12. CD8	Alphanumeric
15	13. CD10	Alphanumeric
	14. CD11c	Alphanumeric
	15. CD13	Alphanumeric
	16. MY4	Alphanumeric
	17. Mo2	Alphanumeric
20	18. CD16	Alphanumeric
	19. CD19	Alphanumeric
	20. CD20	Alphanumeric
	21. CD22	Alphanumeric
	22. CD25	Alphanumeric
25	23. CD33	Alphanumeric
	24. CD34	Alphanumeric
	25. CD41	Alphanumeric
	26. CD56	Alphanumeric
	27. CD57	Alphanumeric
30	28. CD61	Alphanumeric
	29. Kappa	Alphanumeric
	30. Lambda	Alphanumeric
	31. Ia (HLA-DR)	Alphanumeric
	32. Glycophorin A	Alphanumeric
35	33. Cytoplasmic Kappa	Alphanumeric
	34. Cytoplasmic Lambda	Alphanumeric
	35. Tdt	Alphanumeric
	36. Myeloperoxidase Antibody	Alphanumeric
	37. Cytoplasmic CD3	Alphanumeric
40	38. Cytoplasmic CD22	Alphanumeric
	39. Cytoplasmic mu	Alphanumeric
	40. Myeloperoxidase Cytochemistry	Alphanumeric
	41. CD45	Alphanumeric
	42. CD71	Alphanumeric
45	43. Pattern	Numeric
	44. Absolute CD4 Count	Numeric
	45. Size of Abnormal Cells	Alphanumeric

-69-

	46. Percent Abnormal Cells	Numeric
	47. DNA Index	Numeric
	48. Total S-phase	Numeric
	49. Tumor S-phase	Numeric
5	50. Name of Immunophenotyping Panel	Alphanumeric

-70-

## APPENDIX 4: BONE MARROW DATABASE STRUCTURE

	<u>Field Name:</u>	<u>Field Type:</u>
	1. Specimen Identification Number	Alphanumeric
	2. Specimen Date	Paradox Date (long integer)
5	3. Patient Identification Number	Alphanumeric
	4. Specimen Type (aspirate/biopsy/ touch prep)	Alphanumeric
	5. Aspirate Cellularity [%]	Numeric
	6. Differential -- Blasts	Numeric
10	7. Differential -- Promyelocytes	Numeric
	8. Differential -- Myelocytes	Numeric
	9. Differential -- Metamyelocytes	Numeric
	10. Differential -- Neutrophils/Bands	Numeric
	11. Differential -- Monocytes	Numeric
15	12. Differential -- Eosinophils	Numeric
	13. Differential -- Basophils	Numeric
	14. Differential -- Lymphocytes	Numeric
	15. Differential -- Plasma Cells	Numeric
	16. Differential -- Early Erythrocyte Precursors	Numeric
20	17. Differential -- Late Erythrocytes Precursors	Numeric
	18. Other Mononuclear Cells (OMNC)	Numeric
	19. Type of OMNC	Alphanumeric
25	20. Type of Blast Morphology	Alphanumeric
	21. Aspirate Features	Alphanumeric (coded bitmap)
	22. Biopsy Cellularity [%]	Numeric
	23. Biopsy Findings	Alphanumeric (coded bitmap)
	24. Ring Sideroblasts [%]	Numeric
30	25. Non-Specific Esterase Stain [% Positive]	Numeric
	26. Special Stains	Alphanumeric (coded bitmap)
	27. Pattern	Numeric



-71-

# **APPENDIX 5: EXAMPLE CODE FROM PERIPHERAL BLOOD MODULE - PATTERN 7**

```

procedure Pattern7(HWindow:HWnd);
begin
  5   With ACase ^ do
      begin
        StrCat(S,#13#10);
        if (Age<40) or (AbsLymph<15.0) then
          StrCat(S,'The differential diagnosis is a REACTIVE LYMPHOCYTOSIS versus a
10      LYMPHOPROLIFERATIVE DISORDER.'#13#10#13#10);
        if (Age>=40) and (AbsLymph>=15.0) then
          StrCat(S,'The findings suggest a chronic LYMPHOPROLIFERATIVE DISORDER in this
            age group.'#13#10#13#10);
        if (Age<35) and (Ninput[19]='19a') then
15      StrCat(S,'Serologic studies may be helpful to rule-out INFECTIOUS
            MONONUCLEOSIS.'#13#10#13#10);
        if (Age>=40) and (AbsLymph>=15.0) and (CLFLag=FALSE) and (SczFlag=FALSE) and
            (ATLLFLag=FALSE) then
          begin
20      StrCat(S,'The marked absolute lymphocytosis composed of small, "mature-appearing"
            lymphocytes');
            StrCat(S,' is suggestive of CLL.'#13#10#13#10);
            if SmudgeFlag then StrCat(S,'Smudge cells are commonly seen in CLL.'#13#10#13#10);
            end;
25      StrCat(S,'Physical examination for lymphadenopathy and hepatosplenomegaly may be
            helpful.'#13#10#13#10);
            if Ninput[3]<'3d' then
              begin
                if (SpherocyteFlag) and (PolychromFlag) and (SchistoFlag=FALSE) then
30              begin
                StrCat(S,'The presence of spherocytes and polychromasia suggests an immune hemolytic ');
                StrCat(S,'anemia. A COOMB"S test may be helpful.'#13#10#13#10);
                end;
                if (SchistoFlag) and (PolychromFlag) then
35              begin
                StrCat(S,'The presence of schistocytes and polychromasia suggests a microangiopathic ');
                StrCat(S,'hemolytic anemia vs valvular heart disorder. ');
                StrCat(S,'The exact etiology must be determined clinically.'#13#10#13#10);
                end;
40              end;
                if (LGLFlag) and (AbsLymph<15.0) then
                  begin
                    StrCat(S,'The findings suggest a LARGE GRANULAR
                      LYMPHOCYTOSIS.'#13#10#13#10);
45                    StrCat(S,'Clinical history and serologic studies to rule-out rheumatoid arthritis may be
                      helpful.'#13#10#13#10);
                      end;
                      if (Age<40) or (AbsLymph<15.0) then
                        begin
70                        StrCat(S,'If a LYMPHOPROLIFERATIVE DISORDER is suspected,
                            IMMUNOPHENOTYPING by
                            FLOW CYTOMETRY and/or ');
                            StrCat(S,'a BONE MARROW ASPIRATE/BIOPSY may be helpful.'#13#10#13#10);
                            end;
75                        if (Age>=40) and (AbsLymph>=15.0) then
                          begin

```

-72-

StrCat(S, 'The following test may be helpful' #13#10);  
StrCat(S, 'IMMUNOPHENOTYPING by FLOW CYTOMETRY, BONE MARROW  
ASPIRATE/BIOPSY.' #13#10#13#10);

end;

5

end

end;

-73-

# APPENDIX 6: EXAMPLE CODE FROM FLOW CYTOMETRY MODULE FROM PATTERN 8

```

if (8 in PatternMatch) then
begin
5  if (Markers[CD2].Status<>'Not Done') and (Markers[CD3].Status<>'Not Done') and
    (Markers[CD5].Status<>'Not Done') and (Markers[CD7].Status<>'Not Done') then
    begin
      if (Markers[CD2].Status='NEG') or (Markers[CD3].Status='NEG') or
10     (Markers[CD5].Status='NEG')
      or (Markers[CD7].Status='NEG') then AberrantTFlag:=TRUE;
    end;
    if FewMalignant=FALSE then StrCat(IS,'The differential diagnosis includes: '#13#10)
    else StrCat(IS,'The differential diagnosis is suggestive of: '#13#10);
    if (Markers[CD8].Status='POS') and (Markers[CD4].Status='NEG') then
15     StrCat(IS,'T-Cell lymphoproliferative disorders (LGL, T-PLL, ATLL, CTCL), ')
    else
      StrCat(IS,'T-Cell lymphoproliferative disorders (T-PLL, CTCL, ATLL, LGL), ');
      if ((DNAIndex=0.0) or (DNAIndex=1.0)) and (AberrantTFlag=FALSE) then
        StrCat(IS,'Non-Hodgkin's Lymphomas, and Reactive Conditions. ');
20     if ((DNAIndex<>0.0) and (DNAIndex<>1.0)) or (AberrantTFlag) then StrCat(IS,'and Non-
        Hodgkin's
        Lymphomas. ');
        StrCat(IS,'#13#10#13#10);
        if (8 in PatternMatch) and (Markers[CD4].Status='POS') and (Markers[CD8].Status='NEG')
25     and (StrPos(TempString,'Sezary') <> nil) then
          begin
            StrCat(IS,'The helper T-cell phenotype is consistent with the Sezary cells ');
            StrCat(IS,'identified on the peripheral blood smear.'#13#10#13#10);
          end;
30     if AberrantTFlag then
      StrCat(IS,'The findings suggest an aberrant T-phenotype with loss of a pan-T
      marker.'#13#10#13#10);
      if (Markers[CD25].Status='POS') and (Markers[CD4].Status='POS') and
        (StrPos(TempString,'ATLL') = nil)
35     then
        begin
          StrCat(IS,'CD25 positivity is suggestive of ATLL. The HTLV-I status should ');
          StrCat(IS,'be determined to confirm the diagnosis.'#13#10#13#10);
        end;
40     if (Markers[CD25].Status='POS') and (StrPos(TempString,'ATLL') <> nil) then
      begin
        StrCat(IS,'The CD25 positivity is consistent with the presence of circulating "flower" cells. ');
        StrCat(IS,'The HTLV-I status should be determined to confirm the
        diagnosis.'#13#10#13#10);
45     end;
      if (Markers[CD16].Status='POS') or (Markers[CD56].Status='POS')
      or (Markers[CD57].Status='POS') then
        begin
          StrCat(IS,'The positive marker(s) from the group CD16/CD56/CD57 are suggestive of LGL. ');
50     StrCat(IS,'#13#10#13#10);
        end
      end;

```

- 74 -

**APPENDIX 7: EXAMPLE CODE FROM BONE MARROW MODULE**  
**(Myelodysplastic Syndromes)**

```

procedure MDS;
begin
5   With ACase ^ do
    begin
      MDSFlag:=FALSE;
      if (BMBlasts<5) and (RS<=15) and (PBBlasts<=1) and ((PBMono*WBC/100)<=1.0)
      and (Hgb<12) and (AuerFlag=FALSE) and (BMAuerFlag=FALSE) and
10   (FagotFlag=FALSE) and
      (BMFagotFlag=FALSE) then
        begin
          StrCat(S,'Refractory Anemia. ');
          MDSFlag:=TRUE;
15   end;
      if (BMBlasts<5) and (RS<=15) and (PBBlasts<=1) and ((PBMono*WBC/100)<=1.0)
      and (Hgb>=12) and ((WBC<4.0) or (PLT<150)) and (AuerFlag=FALSE) and
      (BMAuerFlag=FALSE)
      and (FagotFlag=FALSE) and (BMFagotFlag=FALSE) then
20   begin
          StrCat(S,'Refractory Cytopenia. ');
          MDSFlag:=TRUE;
          end;
      if (BMBlasts<5) and (RS>15) and (PBBlasts<=1) and ((PBMono*WBC/100)<=1.0)
      and (Hgb<12) and (AuerFlag=FALSE) and (BMAuerFlag=FALSE) and
25   (FagotFlag=FALSE) and (BMFagotFlag=FALSE) then
        begin
          StrCat(S,'Refractory Anemia with Ring Sideroblasts. ');
          MDSFlag:=TRUE;
30   end;
      if (BMBlasts>=5) and (BMBlasts<=20) and (PBBlasts<5) and
      ((PBMono*WBC/100)<=1.0) and (Hgb<12) then
        begin
          StrCat(S,'Refractory Anemia with Excess of Blasts. ');
35   MDSFlag:=TRUE;
          end;
      if (CancerHx[0]<='0') and (BMBlasts<=20) and (PBBlasts<5) and
      (PBMono*WBC/100)>1.0) and (PBBaso*WBC<0.2 and ((PBPro+PBMyel+PBMeta)<=15)
      and (AuerFlag=FALSE)
40   and (BMAuerFlag=FALSE) and (FagotFlag=FALSE) and (BMFagotFlag=FALSE) then
        begin
          if NSE>=20 then StrCat(S,'Chronic myelomonocytic leukemia. ');
          if (NSE<20) and (BMBlasts>=5) and (Hgb<12) then StrCat(S,'RAEB with peripheral
          monocytosis. ');
45   if (NSE<20) and (BMBlasts<5) and (Hgb<12) then StrCat(S,'"Refractory" anemia with
          peripheral
          monocytosis. ');
          if NSE=-1 then
            begin
50   StrCat(S,'#13#10#13#10');
          StrCat(S,'NOTE: ');
          StrCat(S,'A non-specific esterase stain is needed to rule-out CMMoL. ');
          end;
          MDSFlag:=TRUE;
55   end;
    end;
  end;

```

-75-

```
    if ((BMBlasts>20) and (BMBlasts<30)) or (PBBlasts>=5) or (AuerFlag) or (BMAuerFlag)
    and (Hgb<12) then
    begin
    5 StrCat(S,'Refractory Anemia with Excess Blasts in Transformation. ');
      MDSFlag:=TRUE;
    end
    end
    end;
```

-76-

Having thus described the invention, the same will become better understood from the appended claims wherein it is described in a non-limiting manner.

-77-

CLAIMS

1. A method of interpreting hematopathology test results to arrive at a diagnosis of a patient characterized by:

5                   retrieving data from at least one hematopathology test for a patient and storing said hematopathology test data in a database;

                  comparing said hematopathology test data with at least one set of a predetermined number of  
10   hematopathology patterns, each of said hematopathology patterns corresponding to different patient conditions and said patterns arranged in a predetermined hierarchy, to determine which patterns match said hematopathology test data in said  
15   database; and

                  displaying the patterns matched, with recommendations corresponding to the matched patterns for one of no additional testing or one of a plurality of identified additional testing, the  
20   results of which are likely to lead to a diagnosis of the patient's condition.

2. The method according to claim 1 further characterized by:

                  conducting any additional testing  
25   recommended, retrieving the data therefrom and

-78-

storing it in said database with said hematopathology test data;

comparing said additional test data and hematopathology test data with said predetermined  
5 number of patterns, to determine which patterns match said hematopathology test data in said database; and

identifying the highest ranked pattern in said hierarchy to arrive at a differential  
10 diagnosis of the patient's condition.

3. The method according to claim 1 characterized by said hematopathology test being a hematology test in which peripheral blood analysis is conducted.

15 4. The method according to claim 2 characterized by said hematopathology test being a hematology test in which peripheral blood analysis is conducted.

20 5. The method according to claim 4 further characterized by said additional test being a review of a blood smear from the patient, and further including one of entering additional data from the blood smear review and correcting for any erroneous data entered from said hematology test.



-79-

6. The method according to claim 5  
characterized by said differential diagnosis  
including recommendations for additional testing  
including flow cytometry testing of a sample from  
5 the patient, and including:

conducting said flow cytometry testing;

entering the flow cytometry test data  
from said testing into said database with said  
hematology test data;

10 comparing said flow cytometry test data  
and hematology test data with said at least one set  
of a predetermined number of patterns, to determine  
which patterns match said hematology test data and  
flow cytometry test data in said database; and  
15 identifying at least one highest ranked  
pattern in said hierarchy to arrive at a  
differential diagnosis of the patient's condition.

7. The method according to claim 5 further  
characterized by said differential diagnosis  
20 including recommendations for additional testing  
including bone marrow testing of a sample from the  
patient, and further including:

conducting said bone marrow testing;

entering the bone marrow test data from  
25 said testing into said database with said hematology  
test data;

-80-

comparing said bone marrow test data and  
hematology test data with said at least one set of  
a predetermined number of patterns, to determine  
which patterns match said hematology test data and  
5 bone marrow test data in said database; and

identifying the highest ranked pattern  
in said hierarchy to arrive at a differential  
diagnosis of the patient's condition.

8. The method of interpreting  
10 hematopathology test results to arrive at a  
diagnosis of a patient; characterized by:

retrieving data from at least one  
hematopathology test for a patient and storing said  
hematopathology test data in a database;

15 comparing said hematopathology test data  
with at least one set of a predetermined number of  
hematopathology patterns, each of said  
hematopathology patterns corresponding to different  
patient conditions and said patterns arranged in a  
20 predetermined hierarchy to determine which patterns  
match said hematopathology test data in said  
database; and

identifying at least one highest ranked  
pattern in said hierarchy to arrive at a  
25 differential diagnosis of the patient's condition.

-81-

9. The method according to claim 8 further characterized by said hematopathology test being a flow cytometry test, and said comparing step being conducted by comparing said test data from said flow  
5 cytometry test and data from a hematology test previously entered into said database with said patterns.

10. The method according to claim 8 further characterized by said hematopathology test being a  
10 bone marrow test, and said comparing step conducted by comparing said test data from said bone marrow test and data from a hematology test previously entered into said database with said at least one set of patterns.

15 11. The method according to claim 8 further characterized by said hematopathology test being a flow cytometry test, and said comparing step conducted by comparing said test data from said flow  
cytometry test and data from a bone marrow test  
20 previously entered into said database with said at least one set of patterns.

12. The method according to claim 8 further characterized by said hematopathology test being a hematology test, and said comparing step conducted

-82-

by comparing said test data from said hematology test and data from a flow cytometry test previously entered into said database with said at least one set of patterns.

5           13. The method according to claim 8 further characterized by said hematopathology test being a hematology test, and said comparing step conducted by comparing said test data from said hematology test and data from a bone marrow test previously  
10 entered into said database with said patterns.

          14. The method according to claim 8 further characterized by said hematopathology test being a bone marrow test, and said comparing step conducted by comparing said test data from said bone marrow  
15 test and data from a flow cytometry test previously entered into said database with said patterns.

          15. A system for interpreting hematopathology test results to arrive at a diagnosis of a patient, characterized by:  
20           database storage means for storing a database of test results from hematopathology tests for specific patients;

          knowledge base/inference engine means for storing at least one set of a predetermined

-83-

number of patterns corresponding to specific patient conditions, including a normal condition, with said patterns arranged in a predetermined hierarchy, and for storing instructions for comparing said patterns  
5 to said database of test results for determining which patterns are matched;

comparing means for comparing data from said tests for a specific patient with said patterns to determine which patterns are matched and which  
10 pattern is highest in the hierarchy;

input means for inputting said data from test results into said database; and

output means for outputting the results of said comparison by said comparing means, said  
15 results being comprised of at least one of the highest ranked pattern, a differential diagnosis and recommendations for additional testing.

16. The system according to claim 15 further characterized by said input means including a  
20 keyboard, said output means including a display, said database storage means and knowledge base storage means including a "direct access storage device" (DASD), and said comparing means including a central processing unit.

-84-

17. The system according to claim 16 further characterized by said input means further including interface means connected to said comparing means and configured for being connected to a hematology analyzer for having test results therefrom directly  
5 input into said database storage means.

18. The system according to claim 16 further characterized by said input means further including interface means connected to said comparing means  
10 and configured for being connected to a flow cytometry analyzer for having test results therefrom directly input into said database storage means.

19. The system according to claim 15 further characterized by textbook storage means for storing  
15 textbook reference data about hematopathology and connected for allowing a user to make reference to hematopathology reference data when reviewing the results of any comparison conducted with said system.

20. The system according to claim 16 further characterized by textbook storage means for storing  
20 textbook reference data about hematopathology and connected for allowing a user to make reference to hematopathology reference data when reviewing the

-85-

results of any comparison conducted with said system and said textbook storage means including a "direct access storage device" (DASD).

1/20

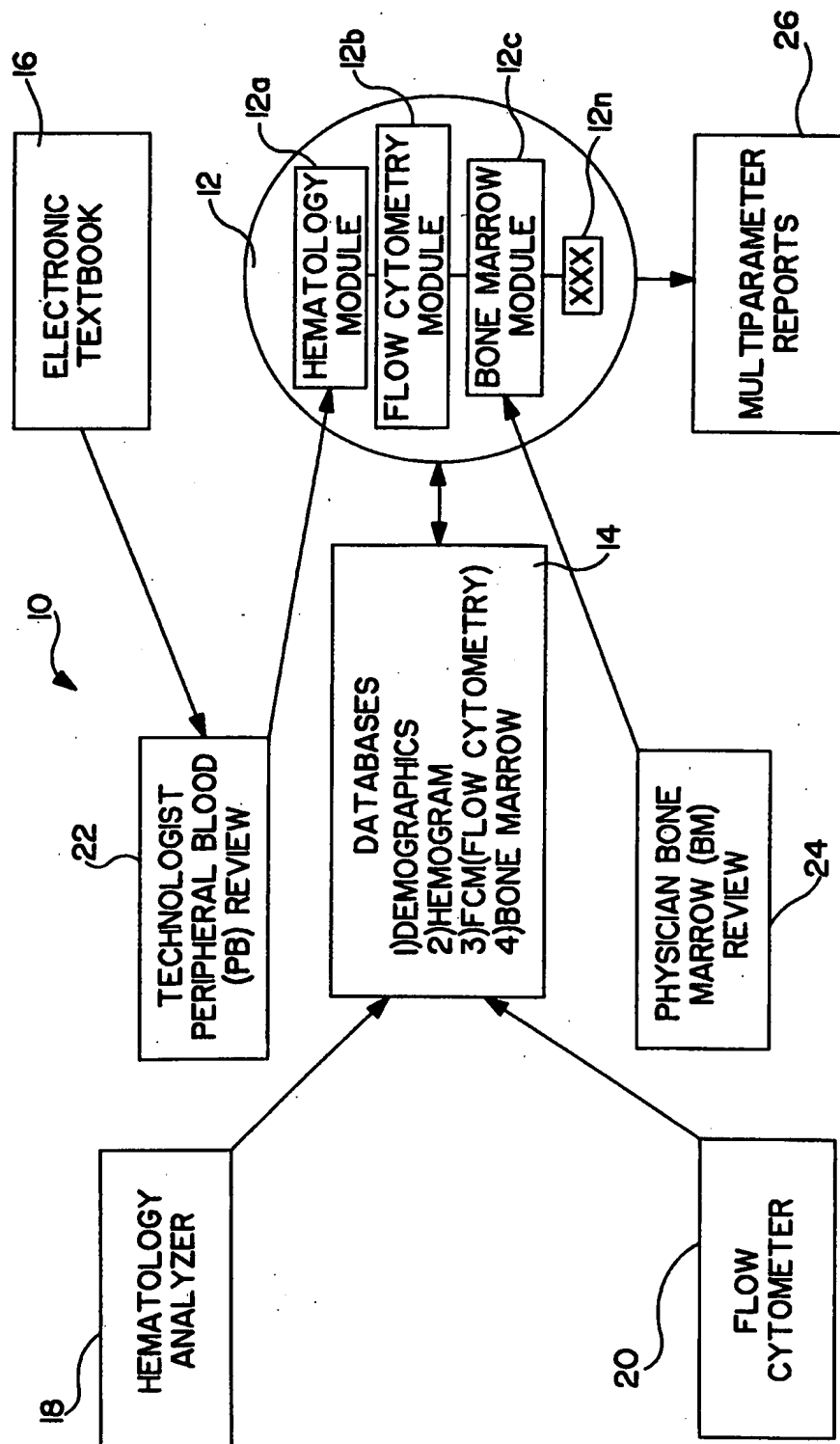


FIG. 1



2/20

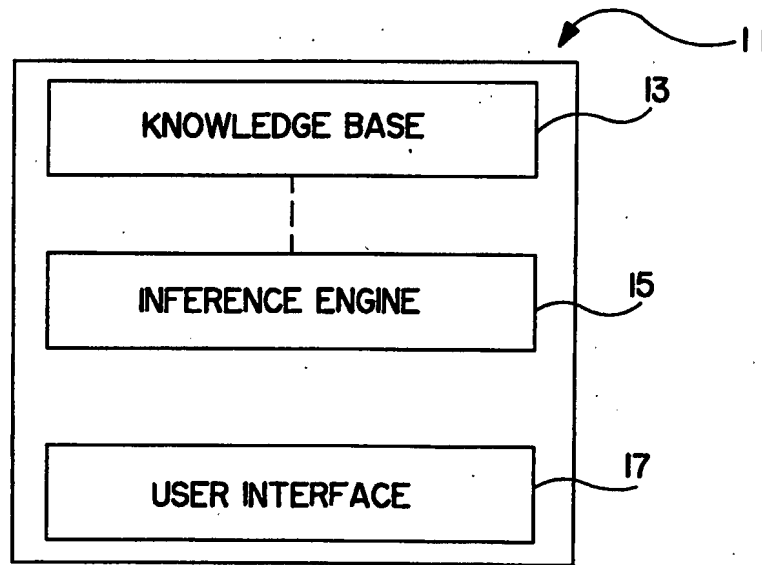


FIG. 2

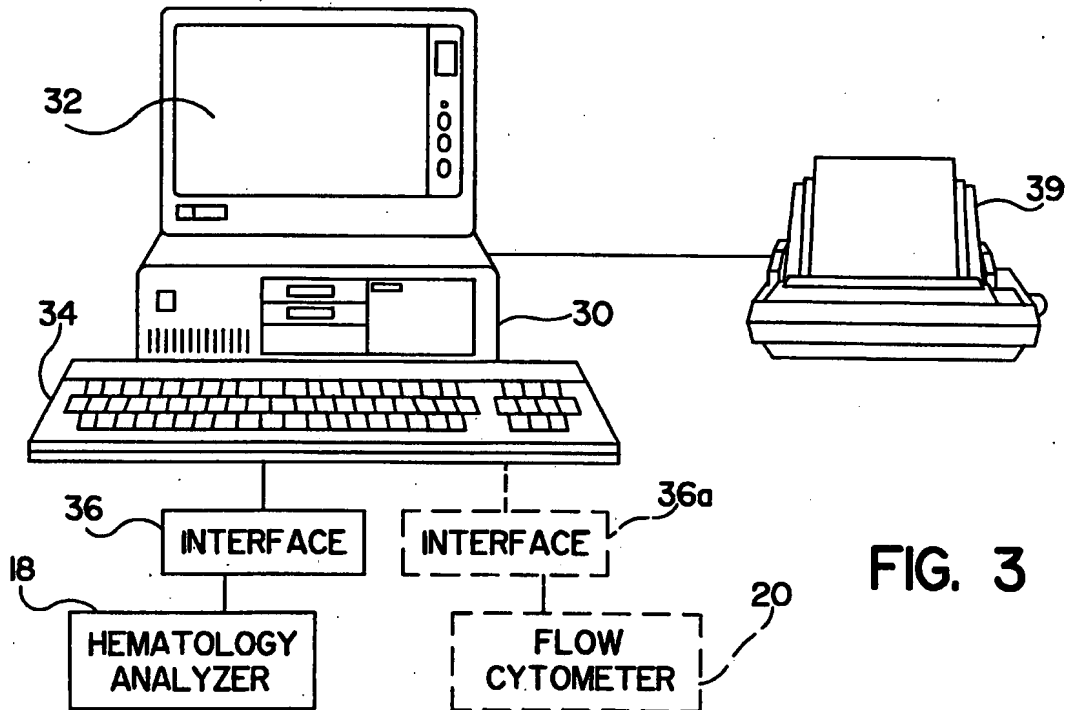


FIG. 3

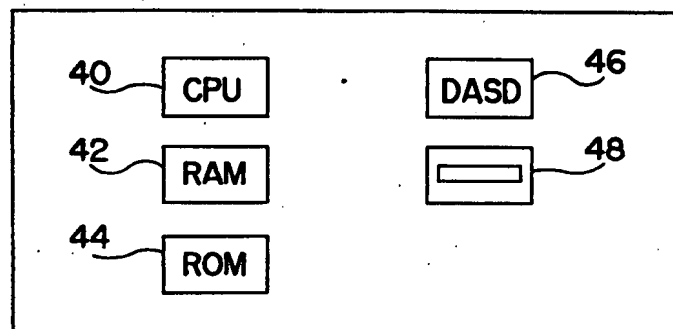


FIG. 4

3/20

AGE: 87  
SEX: F

ID NUMBER: 09990007  
SEQUENCE #: 0999  
DATE OF BIRTH: 3/3/1905  
LOCATION: HHM 752  
PHYSICIAN: Dr. F.C.

COMMENTS:

### Suspect Flags

- ☐ BLASTS
- ☐ IMMATURE GRANULOCYTES/BANDS
- ☐ VARIANT LYMPHS
- ☒ NRBCs
- ☒ DIMORPHIC RBC
- ☒ MICRO RBCs/RBC FRAGMENTS
- ☐ RBC AGGLUTINATION
- ☐ PLATELET CLUMPS

ENTER MANUAL COUNTS

PROCEED

CHOOSE NEW CASE

NAME: C.M.  
DATE/TIME: 3/26/1992 12:12:12  
CASS/POS: 000003  
USER FIELD 1:  
USER FIELD 2:  
USER FIELD 3:

RBC: 2.55  
Hgb: 8.4  
Hct: 26.7  
MCV: 104.6  
MCH: 32.8  
MCHC: 31.4  
RDW: 38.2  
PLT: 668  
WBC: 7.4  
Neut: 5.2  
Lymph: 1.4  
Mono: 0.7  
Eos: 0.1  
Baso: 0.0  
WBC POPULATION: NORMAL  
RBC POPULATION: ABNORMAL  
PLT POPULATION: ABNORMAL

FIG. 5

4/20

Professor Petrushka's Preliminary Comments

Case Identification... 09990007 C.M.

In this 87 year old female the preliminary hemogram findings are consistent with a MACROCYTIC PATTERN with Moderate Normochromic Macrocytic Anemia Marked Anisocytosis Mild Thrombocytosis

Scan the peripheral blood smear (20X objective). Pay attention to RBC and WBC morphology (50X objective). If abnormal cells are seen on scan, a manual WBC differential is recommended.

RBC

☐ Agglutination
☐ Rouleaux
☒ NRBCs

WBC

☐ Accept automated differential?

☐ Yes
☒ No

PLT

☐ Platelet Clumping

Morphologic review necessary?

☒ RBC
☒ WBC
☐ PLT

A manual differential is recommended when many NRBCs are present

Proceed

FIG. 6

5/20

Professor Petrushka's Preliminary Comments

Case identification... 09990007 C.M.

In this 87 year old female the preliminary hemogram findings are consistent with a MACROCYTIC PATTERN with Moderate Normochromic Macrocytic Anemia Marked Anisocytosis Mild Thrombocytosis

Scan the peripheral blood smear (20X objective). Pay attention to RBC and WBC morphology (50X objective). If abnormal cells are seen on scan, a manual WBC differential is recommended.

RBC ☐ Agglutination ☐ Rouleaux ☐ NRBCs

WBC ☐ Accept automated differential? ☒ Yes ☐ No

PLT ☐ Platelet Clumping

Morphologic review necessary? ☒ RBC ☐ WBC ☐ PLT

Proceed

FIG. 7

6/20

NAME:	D.H.	ID NUMBER:	07603942	AGE:	73
DATE/TIME:	4/22/1990 09:41:08	SEQUENCE #:	0417	SEX:	F
CASS/POS:	001108	DATE OF BIRTH:	3/12/1917		
USER FIELD 1:		LOCATION:	HIM 223		
USER FIELD 2:		PHYSICIAN:	Dr. R.C.		
USER FIELD 3:		COMMENTS:			

RBC:	2.74	<b>Suspect Flags</b>
Hgb:	9.0	<input checked="" type="checkbox"/> BLASTS
Hct:	26.9	<input checked="" type="checkbox"/> IMMATURE GRANULOCYTES/BANDS I
MCV:	98.1	<input checked="" type="checkbox"/> VARIANT LYMPHS
MCH:	32.7	<input type="checkbox"/> NRBCs
MCHC:	33.4	<input type="checkbox"/> DIMORPHIC RBC
RDW:	14.0	<input type="checkbox"/> MICRO RBCs/RBC FRAGMENTS
PLT:	262	<input type="checkbox"/> RBC AGGLUTINATION
WBC:	25.5	<input type="checkbox"/> PLATELET CLUMPS
Neut:	0.7	
Lymph:	23.5	
Mono:	11	
Eos:	0.1	
Baso:	0.1	
WBC POPULATION:	ABNORMAL	<div>ENTER MANUAL COUNTS</div>
RBC POPULATION:	ABNORMAL	<div>PROCEED</div>
PLT POPULATION:	NORMAL	<div>CHOOSE NEW CASE</div>

FIG. 8

7/20

## Professor Petrushka's Preliminary Comments

Case identification... 07603942 D.H.

In this 73 year old female the preliminary hemogram findings are consistent with an ABNORMAL MONONUCLEAR CELL PATTERN with

Moderate Normochromic Normocytic Anemia

Mild Leukocytosis

Neutropenia (without Leukopenia)

Severe Lymphocytosis

Monocytosis

Scan the peripheral blood smear (20X objective).

Pay attention to RBC and WBC morphology (50X objective).

A manual WBC differential is recommended.

If no blasts or variant lymphs are seen on scan, the manual differential may not be necessary.

RBC

☐ Agglutination☐ Rouleaux☐ NRBCs

WBC

☐ Accept automated differential?☐ Yes☒ No

PLT

☐ Platelet Clumping

Morphologic review necessary?

☒ RBC☒ WBC☐ PLT

Proceed

FIG. 9

8/20

Spherocytes	<input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+	Polychromasia	<input type="radio"/> 0 <input checked="" type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+
Target Cells	<input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+	Hypochromia	<input checked="" type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+
Teardrop Cells	<input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+	Microcytes	<input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+
Sickle Cells	<input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+	Macrocytes	<input type="radio"/> 0 <input checked="" type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+
Schistocytes	<input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+	Basophilic Stippling	<input checked="" type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+
Ellip/Ovalocytes	<input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+	Howell Jolly	<input checked="" type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+
Stomatocytes	<input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+	Dimorphism	<input type="checkbox"/> Present
Acanthocytes	<input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+	Double population	<input type="checkbox"/> Present
Echinocytes	<input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+	Intracellular organisms	<input type="checkbox"/> Present
Bite Cells	<input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+	Extracellular organisms	<input type="checkbox"/> Present
WBC or PLT Morphology/Differential Required? <input checked="" type="radio"/> Yes <input type="radio"/> No		Proceed	

FIG. 10

9/20

Neutrophils	3				<input type="checkbox"/> Present
Bands					<input type="checkbox"/> Present
Lymphocytes	11				<input type="checkbox"/> Present
Monocytes					<input type="checkbox"/> Present
Eosinophils					<input type="checkbox"/> Present
Basophils					<input type="checkbox"/> Present
NRBCs					<input type="checkbox"/> Present
Blasts	86				<input type="checkbox"/> Present
Promyelocytes					<input type="checkbox"/> Present
Myelocytes					<input type="checkbox"/> Present
Metamyelocytes					<input type="checkbox"/> Present
Other MNCs					<input type="checkbox"/> Present

Hyposegmentation	<input type="checkbox"/> Present
Hypersegmentation	<input type="checkbox"/> Present
Hypogranulation	<input type="checkbox"/> Present
Toxic Granulation	<input type="checkbox"/> Present
Auer Rods	<input type="checkbox"/> Present
WBC inclusions	<input type="checkbox"/> Present
Giant Platelets	<input type="checkbox"/> Present
Micromegakaryocytes	<input type="checkbox"/> Present

Manual Platelet

<input type="checkbox"/> Variant lymphs	<input type="checkbox"/> Plasmacytoid lymphs	<input type="checkbox"/> Large granular lymphs	<input type="checkbox"/> Smudge Cells
---	--	--	---------------------------------------

Review Morphology	<input checked="" type="checkbox"/> WBC	<input type="checkbox"/> PLT
Manual Diff/Count	<input checked="" type="checkbox"/> WBC	<input type="checkbox"/> PLT

**Proceed**

FIG. 11



10/20

# Professor Petrushka's Interpretation and Recommendations

Case identification... 07603942 D.H.

In this 73 year old female the hemogram findings are characteristic of an ABNORMAL MONONUCLEAR CELL PATTERN with Moderate Normochromic Normocytic Anemia Mild Leukocytosis Neutropenia (without Leukopenia)

The differential diagnosis includes: ACUTE LEUKEMIA.

The presence of blasts with Auer rods is highly suggestive of an ACUTE MYELOID LEUKEMIA.

The following tests may be helpful: BONE MARROW ASPIRATE/BIOPSY with CYTOGENETIC analysis. CYTOCHEMICAL STAINS (Myeloperoxidase, NSE) on blood or bone marrow smears. IMMUNOPHENOTYPING (by FLOW CYTOMETRY) on blood or bone marrow.

Next Case

Choose New Case

FIG. 12

11/20

Flow Cytometry Data Entry									
Case #	<input type="text"/>	Date (mm/dd/yy)	<input type="text"/>	Patient ID#	<input type="text"/>				
Name	<input type="text"/>			Sex	<input type="text"/>	DOB (mm/dd/yy)	<input type="text"/>		
Specimen	<input type="radio"/> PB	<input type="radio"/> BM	<input type="radio"/> LN	<input type="radio"/> FL	<input type="radio"/> SP	<input type="radio"/> EN			
CD1a	<input type="text"/>	CD19	<input type="text"/>	Gly A	<input type="text"/>	% Blasts	<input type="text"/>		
CD2	<input type="text"/>	CD20	<input type="text"/>	cKappa	<input type="text"/>	<div>Abn Mononuclear Cells</div> <div>Burkitt's Cells</div> <div>Cerebriiform Cells</div> <div>Flower Cells</div> <div>Hairy Cells</div> <div>Lymphoma Cells</div> <div>None</div> <div>Plasma Cells</div> <div>Plasmablasts</div> <div>Prolymphocytes</div>			
CD3	<input type="text"/>	CD22	<input type="text"/>	cLambda	<input type="text"/>				
CD4	<input type="text"/>	CD25	<input type="text"/>	Tdt	<input type="text"/>				
CD5	<input type="text"/>	CD33	<input type="text"/>	MPO Ab	<input type="text"/>				
CD7	<input type="text"/>	CD34	<input type="text"/>	cCD3	<input type="text"/>				
CD8	<input type="text"/>	CD41	<input type="text"/>	cCD22	<input type="text"/>				
CD10	<input type="text"/>	CD56	<input type="text"/>	c mu	<input type="text"/>				
CD11c	<input type="text"/>	CD57	<input type="text"/>	MPO Cyto	<input type="text"/>				
CD13	<input type="text"/>	CD61	<input type="text"/>	CD45	<input type="text"/>				
MY4	<input type="text"/>	Kappa	<input type="text"/>	CD71	<input type="text"/>				
Mo2	<input type="text"/>	Lambda	<input type="text"/>						
CD16	<input type="text"/>	Ig	<input type="text"/>						
DNA Index <input type="text"/>		Total S% <input type="text"/>		% Abn Cells <input type="text"/>		Cell Size <input type="text"/>			
		Tumor S% <input type="text"/>							
				OK		Cancel			

FIG. 13

12/20

Name:	F.B.	ID Number:	SMO272361	Additional Information
Sex:	M	Date of Birth:	1/1/1924	
Specimen: LN		Specimen Date:	2/27/1992	
Specimen Number: RB11250				

CD1a:	-	CD61:	-	<p><b>Professor Fidello's Interpretation:</b></p> <p>The immunophenotype is characteristic of a MATURE B-CELL PATTERN MONOCLONAL for LAMBDA.</p> <p>In a lymph node, the differential diagnosis for this pattern includes: B-cell non-Hodgkin's Lymphomas.</p> <p>CD5 positivity and weak surface immunoglobulin staining are highly suggestive of small lymphocytic lymphoma/CLL.</p> <p>The small cell size and low S-phase are consistent with the above diagnosis.</p> <p>The results should be correlated with the lymph node biopsy, clinical features, and appropriate laboratory studies.</p>
CD2:	-	Kappa:	+	
CD3:	-	Lambda:	+	
CD4:	-	Ia(HLA-DR):	+	
CD5:	+	Glycophorin A:		
CD7:	-	cKappa:		
CD8:	-	cLambda:		
CD10:	-	Tdt:		
CD11c:	-	MPO Ab:		
CD13:	-	cCD3:		
MY4:	-	cCD22:		
Mo2:	-	c mu:		
CD16:	++	MPO cyto:	++	
CD19:	++	CD45:	+	
CD20:	+	CD71:	Small	
CD22:	+	Size Abn Cells:	61	
CD25:	-	%Abn Cells:	1.0	
CD33:	-	DNA Index:	1.0	
CD34:	-	Total S%:		
CD41:	-	Tumor S%:		
CD56:	-			
CD57:	-			

Next Case

FIG. 14

13/20

Name:	N.B.	ID Number:	HHMO97	Additional Information L3
Sex:	M	Date of Birth:	4/10/1976	
Specimen:	PB	Specimen Date:	3/15/1990	
Specimen Number:	EX004			

CD1a:	1	CD61:	83	<p><b>Professor Fidelio's Interpretation:</b></p> <p>The immunophenotype is characteristic of a MATURE B-CELL PATTERN MONOCLONAL for KAPPA.</p> <p>In the peripheral blood, the differential diagnosis for this pattern includes: B-ALL(L3) and B-CELL lymphoproliferative disorders.</p> <p>The combination of a mature B-cell phenotype and blasts with Burkitt's cell morphology is consistent with a diagnosis of B-ALL(L3)</p> <p>The results should be confirmed with a bone marrow aspirate/biopsy and correlated with the clinical features and appropriate laboratory studies.</p>
CD2:	7	Kappa:	3	
CD3:	8	Lambda:	83	
CD4:	4	Ia(HLA-DR):		
CD5:	5	Glycophorin A:		
CD7:	3	cKappa:		
CD8:	4	cLambda:		
CD10:		Tdt:		
CD11c:		MPO Ab:		
CD13:	2	cCD3:		
MY4:	4	cCD22:		
Mo2:		c mu:		
CD16:		MPO cyto:		
CD19:	70	CD45:		
CD20:	68	CD71:		
CD22:		Size Abn Cells:		
CD25:		%Abn Cells:		
CD33:	1	DNA Index:		
CD34:	1	Total S%:		
CD41:	1	Tumor S%:		
CD56:				
CD57:				

Next Case

FIG. 15

14/20

Name: B.D.  
 Sex: F  
 Physician: ID Number: HHM162 Date: 3/9/93  
 Case #: BOOI Age: 83  
 History:

Specimen ☒ Aspirate ☐ Biopsy ☐ Touch Prep

## Peripheral Blood Findings

Hgb:	10.0	Neut:	13%	Blasts:	70%
MCV:	95.0	Bands:		Pro:	
MCH:	33.0	Lymph:	15%	Myel:	1%
RDW:	19.6	Mono:		Meta:	1%
PLT:	61	Eos:		MNCs:	
WBC:	159.0	Baso:		Morph:	

## Recommendations

Peripheral Blood(12/13/1988):ABNORMAL MONONUCLEAR CELL PATTERN.  
 In the presence of circulating blasts, an ACUTE LEUKEMIA must be excluded.  
 If the bone marrow demonstrates an ACUTE LEUKEMIA, the following studies are  
 recommended: CYTOCHEMISTRIES, IMMUNOPHENOTYPING, CYTOGENETICS. If  
 the diagnosis is ALL, DNA ANALYSIS may be helpful.

Approve CBC Results?

☒ Yes ☐ No

Choose New Case

Proceed

FIG. 16

15/20

**Low Power Evaluation**

Spicules ☒ Present ☐ Absent

Megakaryocytes ☐ Normal ☒ Decreased ☐ Increased

Cellularity(%)           100%

Foreign Cells

(50X Oil Objective)  Differential

Blasts: 65 Monocytes: Early erythroid: 16  
 Promyelocytes: Eosinophils: Late erythroid: 16  
 Myelocytes: Basophils: MNC (%):  
 Metamyelocytes: Lymphocytes: 18 MNC Type:  
 Neutrophils/Bands: 1 Plasma Cells: ☐ Organisms

Qualitative Abnormalities Monocyte-like morphology

☐ Erythroid Karyorrhexis ☐ Hypogranulation ☐ Abnormal Eosinophils

N/C Asynchrony  Giant Precursors

Cytoplasmic Vacuoles

Special Stains  Proceed

FIG. 17

16/20

<b>Aspirate</b>	
Iron Storage _____ <input type="radio"/> Normal <input type="radio"/> Decreased <input type="radio"/> Increased	
Iron Incorporation _____ <input type="radio"/> Normal <input type="radio"/> Decreased <input type="radio"/> Increased	
Ring Sideroblasts (%)	Not Done
<b>Stains for Blasts and Abnormal Mononuclear Cells</b>	
Peroxidase _____ <input type="radio"/> Negative <input checked="" type="radio"/> Positive	NSE <div style="border: 1px solid black; height: 40px; width: 100%; position: relative;"> <span style="position: absolute; left: -10px; top: 50%; transform: translateY(-50%);">↔</span> <span style="position: absolute; right: -10px; top: 50%; transform: translateY(-50%);">→</span> </div>
<input type="checkbox"/> PAS Positive <input type="checkbox"/> TRAP Positive	
<b>Biopsy</b>	
Reticulin _____ <input type="radio"/> 0 <input type="radio"/> 1+ <input type="radio"/> 2+ <input type="radio"/> 3+	
Positive for Organisms	Proceed

Fig. 18

17/20

**Professor Belmonte's Interpretation and Recommendations**

Peripheral Blood: Hgb 10.0, MCV 95.0, MCH 33.0, PLT 61, WBC 159.0  
 Manual Differential: 13% Neutrophils, 15% Lymphocytes, 70% Blasts, 1%  
 Myelocytes, 1% Metamyelocytes.

In this 83 year old female, the bone marrow pattern is **MONONUCLEAR CELL INFILTRATION**.

The bone marrow aspirate is hypercellular. Megakaryocytes are decreased.

Bone Marrow Differential: 65% Blasts, 1% Neutrophils/Bands, 18%  
 Lymphocytes, 16% late Erythroid precursors.

The blasts are peroxidase positive and have abundant blue, vacuolated cytoplasm  
 consistent with monoblasts. Positive NSE staining was seen in 90% of the cells.

Immunology: The blasts in the peripheral blood demonstrated a **MYELOID**  
 phenotype.

**INTERPRETATION:**

Bone Marrow, aspirate: **Acute Myeloid Leukemia, FAB M5a.**

[Next Case](#)
[Choose New Case](#)
**FIG. 19**



18/20

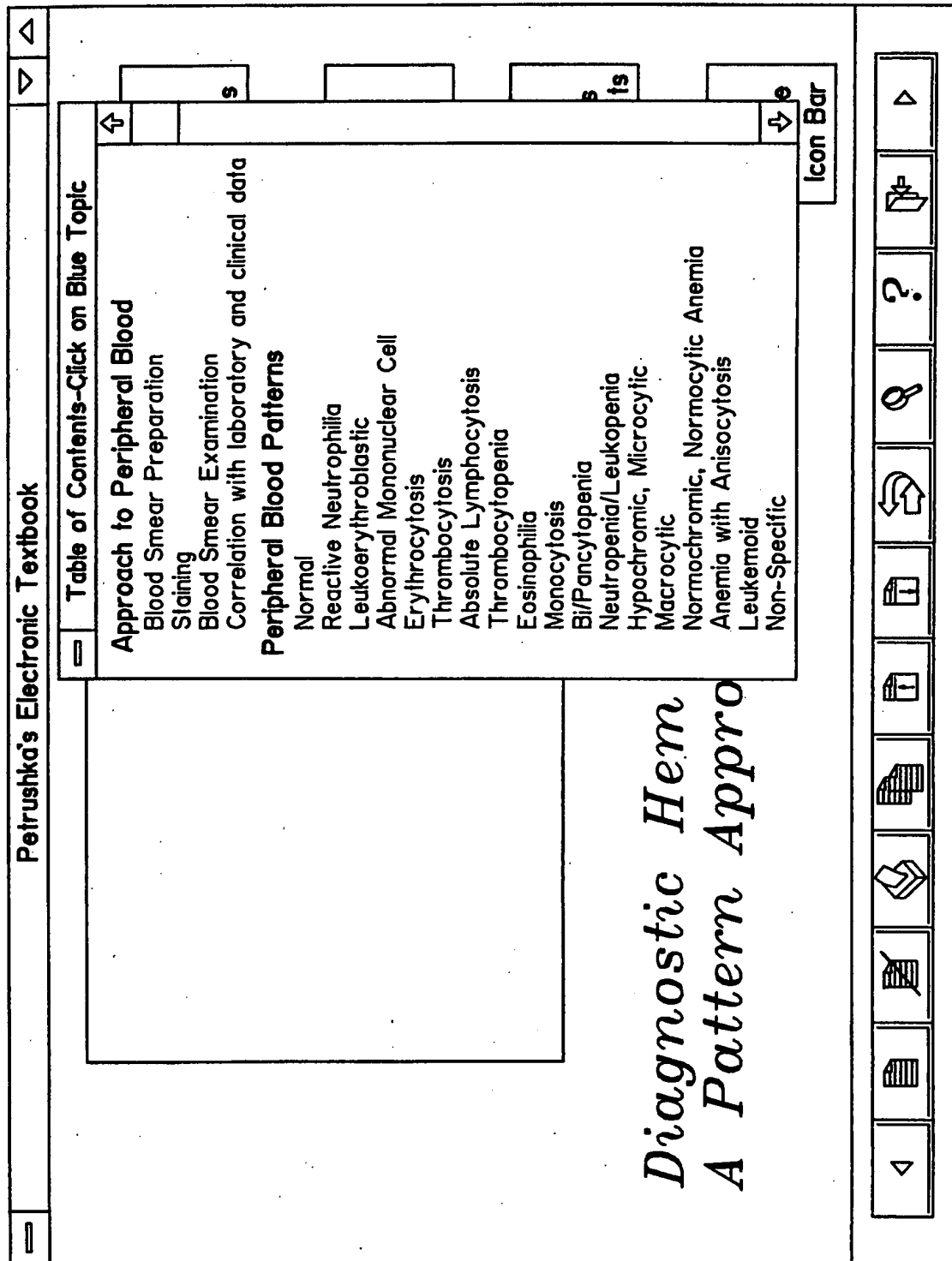
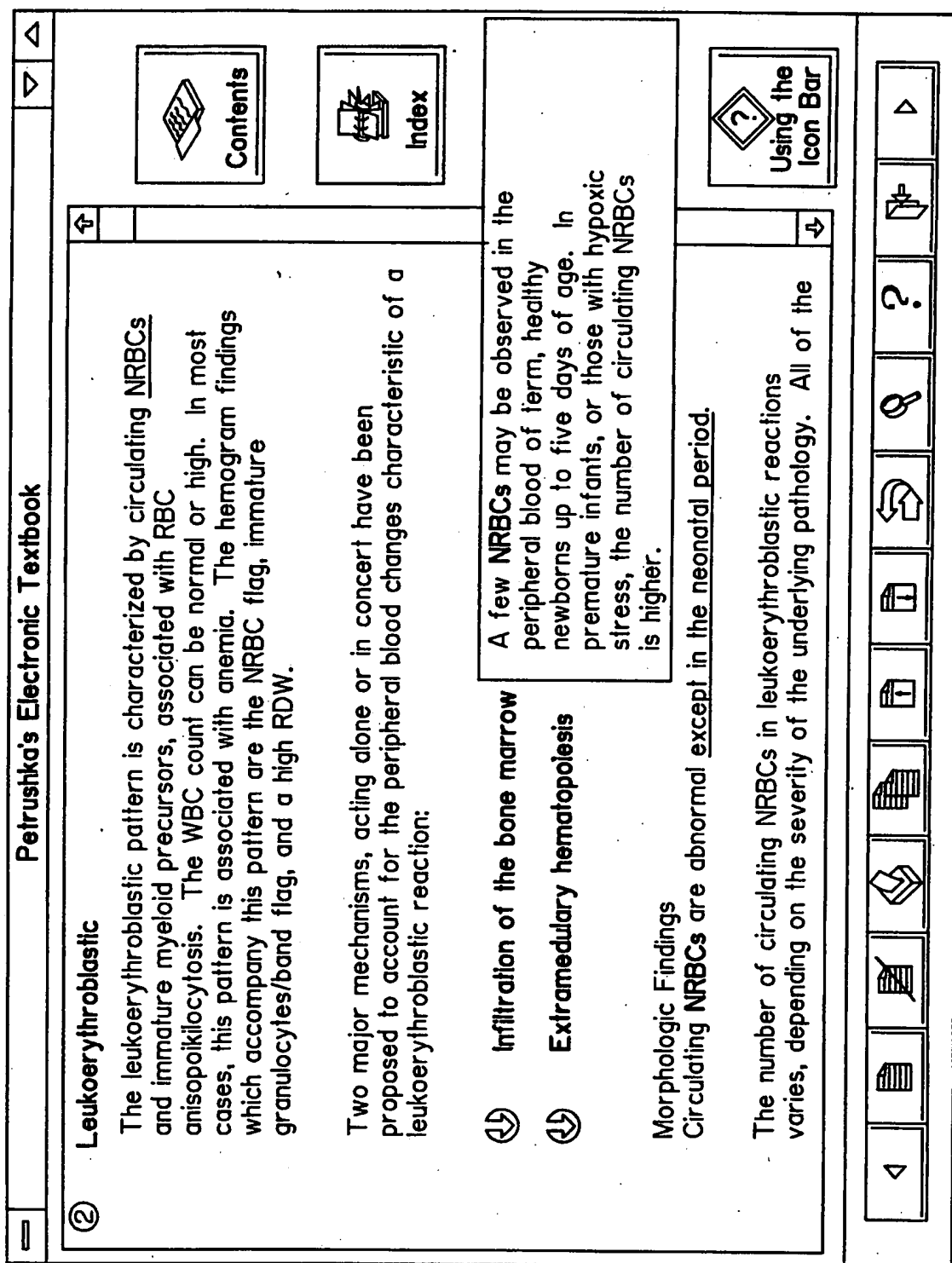


FIG. 20


SUBSTITUTE SHEET (RULE 26)



**FIG. 21**

20/20

FIG. 22

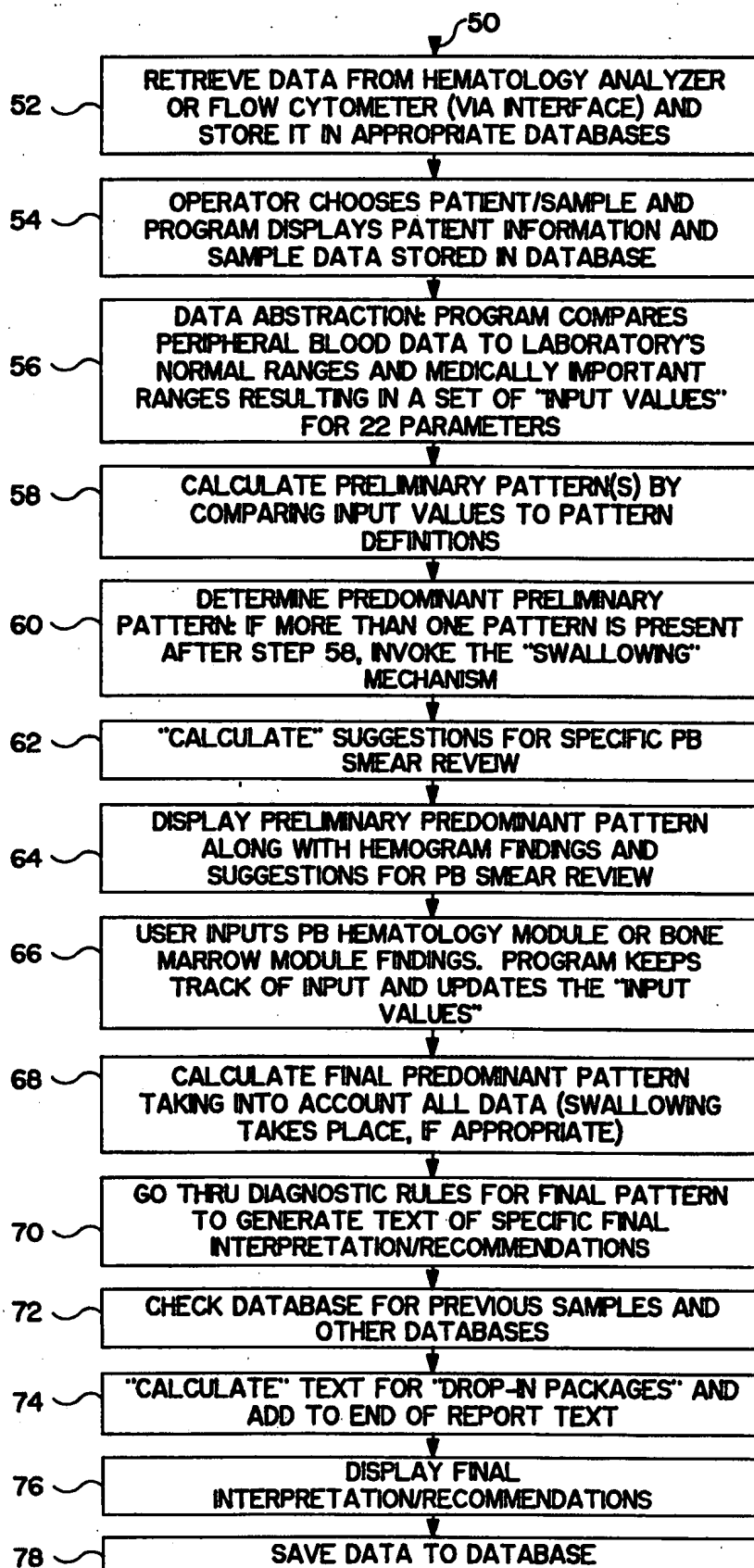

## INTERNATIONAL SEARCH REPORT

International application No.  
PCT/US94/10001**A. CLASSIFICATION OF SUBJECT MATTER**

IPC(6) : Please See Extra Sheet.

US CL : 364/413.07, 413.08, 413.02; 128/630; 436/63, 17, 172; 382/6

According to International Patent Classification (IPC) or to both national classification and IPC

**B. FIELDS SEARCHED**

Minimum documentation searched (classification system followed by classification symbols)

U.S. : 364/413.07, 413.08, 413.02; 128/630; 436/63, 17, 172; 382/6

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practicable, search terms used)

APS

**C. DOCUMENTS CONSIDERED TO BE RELEVANT**

Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	MEDINFO 92, ELSEVIER SCIENCE PUBLISHERS B.V. (NORTH-HOLLAND), ENTITLED "A DECISION SUPPORT SYSTEM FOR DIAGNOSTIC CONSULTATION IN LABORATORY HEMATOLOGY", D. T. NGUYEN ET AL, PAGES 591-595, 05 SEPTEMBER 1992	1 AND 3
Y		2, 4-5, 8, 15-17, AND 19-20
Y	"EVALUATION OF A KNOWLEDGE-BASED SYSTEM FOR INTERPRETING FLOW CYTOMETRIC IMMUNOPHENOTYPING DATA" L. W. DIAMOND ET AL, 1993, pp. 124-128.	6, 9, 12, AND 18
Y	"COMPUTER-ASSISTED BONE MARROW INTERPRETATION: A PATTERN APPROACH", PROCEEDINGS OF MIE 93, APRIL 1993, NGUYEN ET AL, PAGES 119-123.	7, 10, 11, 13, AND 14

☐ Further documents are listed in the continuation of Box C.☐ See patent family annex.

* Special categories of cited documents:	"T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
"A" document defining the general state of the art which is not considered to be part of particular relevance	"X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
"E" earlier document published on or after the international filing date	"Y" document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art
"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)	"Z" document member of the same patent family
"O" document referring to an oral disclosure, use, exhibition or other means	
"P" document published prior to the international filing date but later than the priority date claimed	

Date of the actual completion of the international search

12 NOVEMBER 1994

Date of mailing of the international search report

JAN 23 1995

Name and mailing address of the ISA/US  
Commissioner of Patents and Trademarks  
Box PCT  
Washington, D.C. 20231

Facsimile No. N/A

Authorized officer

XUONG M. CHUNG-TRANS

Telephone No. (703) 305-9711

# INTERNATIONAL SEARCH REPORT

International application No.  
PCT/US94/10001

## A. CLASSIFICATION OF SUBJECT MATTER:

IPC (6):

G06F 15/00; A61B 5/00; G01N 33/48, 31/00, 33/00; G06K 9/00